

i2B2 Webclient Help Documentation

Version 1.8.0

Updated 8/2023

Table of Contents

1. Quick Start	4
1.1 Basic Query	4
[How to]	4
[Example]	4
1.2 Adding an AND term	5
[Example]	6
1.3 Adding an OR term.....	8
[How to]	8
[Example]	8
2. Webclient Layout.....	10
Find Patients (Top Right)	10
Terms (Top Left)	10
Query Status (Bottom Right)	10
Queries (Bottom Left).....	10
Workplace (Middle Left).....	11
3. Query Tool	12
3.1 Layout	12
3.2 Creating a Query	12
3.2.1 i2b2 Items for Querying.....	12
3.2.1 Operators and Order of Precedence	13
3.2.2 Constraints	15
3.2.3 Query Timing.....	23
3.2.4 Reloading a Previous Query	27
3.3 Running a Query	28
3.3.1 Run Query	28
3.3.2 Query Result Options.....	30
3.3.3. Query Status	31
4. Navigate Terms	33

Navigate Terms Tab	33
Find Mode	35
A. The Search by Names mode is used to find all query items containing particular text.....	35
B. The Search by Codes tab is used to find coded variables.....	37
Options for Finding Terms	38
5. Query Status	39
Show Query Status	40
Graph Results	40
Query Report	42
6. Previous Queries	43
6.1 'Queries' tab.....	43
6.2 Query Search mode	43
6.3 Layout of a Previous Query Display	44
6.3.1 Previous Query Name.....	44
6.3.2 Results Folder	45
6.4 Navigation	46
6.4.1 Menu.....	46
6.4.1.1 Display.....	46
6.4.1.2 Rename a Previous Query	46
6.4.1.3 Delete / Remove a Previous Query.....	46
6.4.1.4 Refresh the List.....	46
6.4.2 Options.....	46
6.4.2.1 Maximum to Display	47
6.4.2.2 Sort Queries by Name.....	47
6.4.2.3 Sort by Create Date.....	47
6.4.2.4 Auto Refresh.....	48
6.5 Using a Previous Query	48
6.5.1 Draggable Items.....	48
6.5.1.1 Previous Query	48
6.5.1.2 Patient Set.....	48
6.5.1.3 Encounter Set.....	48
6.5.2 Add to 'Find Patients'	49
6.5.2.1 Add Previous Query to 'Find Patients'	49

6.5.2.2 Add Patient Set to the 'Query Tool'	50
6.5.2.3 Add Encounter Set to the 'Query Tool'	50
6.5.3 Add to 'Workplace' panel.....	50
6.5.3.1 Add Previous Query to 'Workplace' Panel	50
6.5.3.2 Add Patient Set to 'Workplace' panel	50
6.5.3.3 Add Encounter Set to 'Workplace' panel	51
7. Workplace	51
7.1 Layout of the 'Workplace' panel	51
7.1.1 Containers	52
7.1.2 Folders	52
7.1.3 Leaves (Items)	52
7.2 Adding a folder and items.....	52
7.2.1 Folders and Items Draggable to Workplace	52
7.2.2 Adding a Folder	53
7.2.3 Adding Items	53
7.2.3.1 Adding Concept, Concept Folder, Modifier, and Modifier Folder.....	53
7.2.3.2 Adding Previous Query, Patient Set, and Encounter Set.....	53
7.3.1. Expand / Collapse a Folder.....	54
7.3.2 Rename a Folder or an Item.....	54
7.3.3 Annotate a Folder or item	55
7.3.4 Delete a Folder or item	55
7.4 Using Items in 'Workplace' Panel.....	55
7.4.1 Items Draggable from Workplace	55

1. Quick Start

1.1 Basic Query

Basic queries are constructed by dragging and dropping terms from the 'Terms' panel into a 'Find Patients' group panel.


[How to]

1. Find a term
 - a. Go to 'Terms' tab
 - b. Click folders to expand the lists of terms
 - c. Select a term
2. Drag and drop to 'Find Patients' group panel
 - a. Drag the selected term into 'Find Patients> Find Patients' group panel
3. Run Query
 - a. Click 'Find Patients> Run Query Options' → 'Run Query'
4. Get Result
 - a. Query running status will be displayed in 'Query Status' section, which is located below 'Find Patients'.
 - b. Query result will be displayed in 'Query Status'
 - i. Number of patient
 - ii. Query name
5. The query and its result will be stored and displayed in 'Queries' tab

[Example]

Find all patients who have
at least one diagnosis of *Acute respiratory infections*.

1. Find *Acute respiratory infections*
 - a. Go to 'Terms' tab
 - b. Click on *Diagnosis* folder to expand → Click on *Respiratory system* folder to expand → Find *Acute respiratory infections*
2. Drag and drop
 - a. Drag the *Acute respiratory infections* term into 'Find Patients > Find Patients' group panel
3. Run Query
 - a. Click 'Find Patients button > Run Query Options' > 'Run Query'
4. Get Results
 - a. See the 'Query Status' section for the query result.
 - b. The query is now stored in the 'Queries' tab

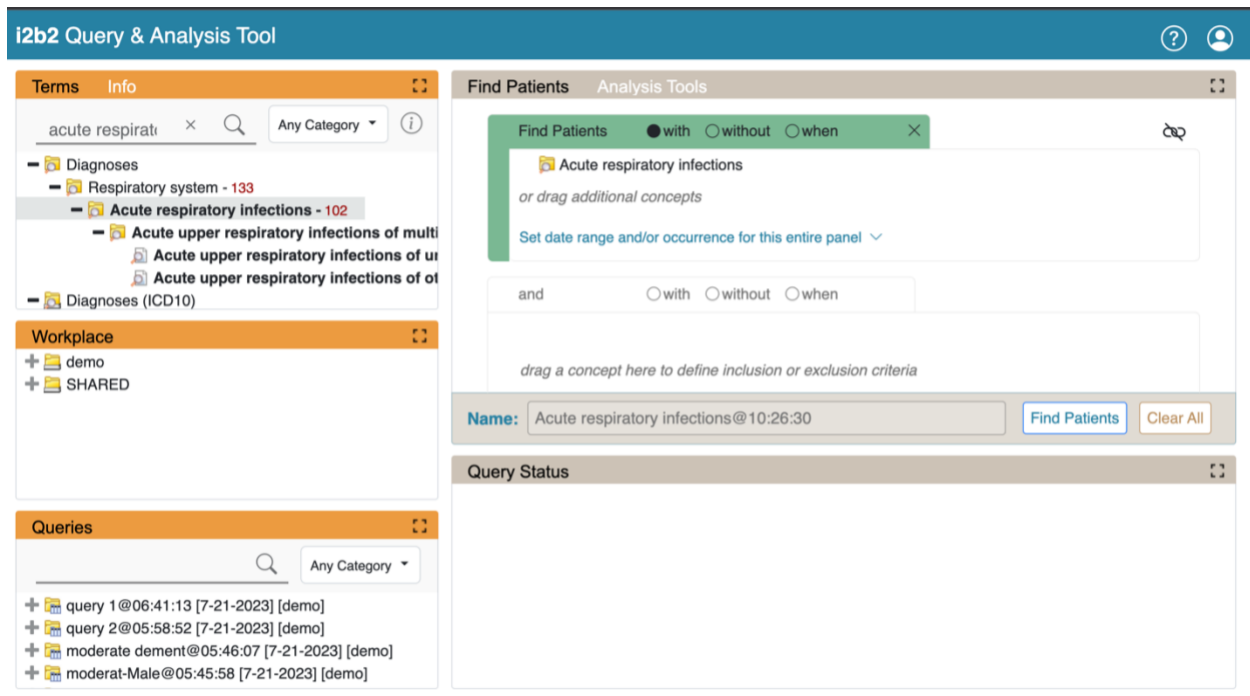

1.2 Adding an AND term

Terms in different groups have a logical AND relationship. When terms are combined with the AND operator, retrieved records must contain all the terms.

[How to]

1. Find a Term
 - a. Go to 'Terms' tab
 - b. Click term folders to narrow down the concept
2. Drag and drop to 'Find Patients' group panel
 - a. Drag the selected term into 'Find Patients > Find Patients' group panel
3. Find another Term
4. Drag and drop to 'and' group panel
 - a. Drag the selected term into 'Find Patients > and' group panel'
5. Run Query
 - a. Click 'Find Patients > Run Query Options' > 'Run Query'
6. Get Result
 - a. See the 'Query Status' tab for the query result
 - b. The query is now stored in the 'Queries' tab

[Example]

Find all patients who have
at least one diagnosis of *Acute respiratory infections*
AND
at least one prescription of *Bronchodilators*.

1. Find *Acute respiratory infections*
 - a. Go to 'Navigate Terms' tab
 - b. Select *Diagnosis* → *Respiratory system* → *Acute respiratory infections*
2. Drag and drop
 - a. Drag the *Acute respiratory infections* term into 'Find Patients > Find Patients' group panel
3. Find *Bronchodilators*
 - a. Go to 'Navigate Terms' panel
 - b. Select *Medications* → *Respiratory agents* → *Bronchodilators*
4. Drag and drop
 - a. Drag the *Acute respiratory infections* term into 'Find Patients> and' group panel'
5. Run Query
 - a. Click 'Find Patients> Run Query Options' → 'Run Query'
6. Get Results
 - a. See the 'Query Status' tab for the query result
 - b. The query is now stored in the 'Queries' tab

Terms Info

Any Category

- Clinical Trials
- Custom Metadata
- Demographics
- Diagnoses
 - Circulatory system - 66
 - Conditions in the perinatal period - 13
 - Congenital anomalies - 23
 - Digestive system - 100
 - Endocrine disorders - 35
 - Events of pregnancy - 14
 - Genitourinary system - 65
 - Hematologic diseases - 27
 - Infectious and parasitic diseases - 86
 - Injury and poisoning - 109
 - Mental Disorders - 81
 - Metabolic and immunity disorders - 78
 - Musculoskeletal and connective tissue - 108
 - Neoplasms - 65
 - Neurologic Disorders - 96
 - Nutritional deficiencies - 3
 - Respiratory system - 133
 - Admit Diagnosis
 - Principal Diagnosis
 - Secondary Diagnosis
 - Acute respiratory infections - 102
 - Chronic obstructive diseases - 133
 - Lung diseases due to external irritants - 2
 - Other diseases - 51
 - Pneumonia and influenza - 50
 - Upper respiratory tract infections - 107
 - Skin diseases - 101
 - Symptoms, signs, and ill-defined conditions - 133
 - E-codes - 55
 - V-codes - 132
- Diagnoses (ICD10)
- Expression Profiles Data
- Laboratory Tests

Find Patients Analysis Tools

Find Patients with without when

Acute respiratory infections

or drag additional concepts

Set date range and/or occurrence for this entire panel

and with without when

drag a concept here to define inclusion or exclusion criteria

Name: Acute respiratory infections@21:39:19

Find Patients Clear All

Query Status

Terms Info

Any Category

- Patient Instructions
- Route
- Rx Sig
- Alternative medicines
- Antihyperlipidemic agents - 7
- Anti-infectives - 78
- Antineoplastics - 3
- Biologicals
- Cardiovascular agents - 19
- Central nervous system agents - 78
- Coagulation modifiers - 11
- Gastrointestinal agents - 43
- Hormones - 49
- Immunologic agents - 3
- Miscellaneous agents - 13
- Nutritional products - 30
- Psychotherapeutic agents - 54
- Respiratory agents - 97
 - Patient Instructions
 - Route
 - Rx Sig
 - Antiasthmatic combinations
 - Antihistamines - 42
 - Antitussives - 8
 - Bronchodilators - 77
 - Decongestants - 4
 - Expectorants - 1
 - Leukotriene modifiers - 15
 - Lung surfactants
 - Miscellaneous respiratory agents
 - Respiratory inhalant products - 72
 - Upper respiratory combinations - 16
 - Topical agents - 66
- Procedures
- Providers
- Reports
- Visit Details

Find Patients Analysis Tools

Find Patients with without when

Acute respiratory infections

or drag additional concepts

Set date range and/or occurrence for this entire panel

and with without when

Bronchodilators

or drag additional concepts

Set date range and/or occurrence for this entire panel

and with without when

drag a concept here to define inclusion or exclusion criteria

Name: Acute respiratory infections-Bronchodilators@21:40:36

Find Patients Clear All

Query Status

1.3 Adding an OR term

Terms in the same group have a logical **OR** relationship. When terms are combined with the **OR** operator, retrieved records will contain either or both terms.

[How to]

1. Find a Term
 - a. Go to 'Terms' tab
 - b. Click term folders to narrow down the concept
 - c. Select a term
2. Drag and drop to 'Find Patients' group panel
 - a. Drag the selected term into 'Find Patients> Find Patients' group panel
3. Find another Term
4. Drag and drop to 'Find Patients' group panel
 - a. Drag the selected term into 'Find Patients> Find Patients' group panel, again.
5. Run Query
 - a. Click 'Find Patients> Run Query Options' → 'Run Query'
6. Get Result
 - a. See the 'Query Status' tab for the query result
 - b. The query is now stored in the 'Queries' tab

[Example]

Find all patients who have

at least one diagnosis of *Acute respiratory infections* OR at least one diagnosis of *Asthma*


AND

at least one prescription of *Bronchodilators*.

1. Find *Acute respiratory infections*
 - a. Go to 'Terms' tab
 - b. Select *Diagnosis* → *Respiratory system* → *Acute respiratory infections*
2. Drag and drop
 - a. Drag the *Acute respiratory infections* term into 'Find Patients> Find Patients' group panel
3. Find *Asthma*
 - a. Go to 'Terms' tab
 - b. Select *Diagnosis* → *Respiratory system* → *Chronic obstructive diseases* → *Asthma*
4. Drag and drop
 - a. Drag the *Acute respiratory infections* term into 'Find Patients> Find Patients' group panel, again
5. Find *Bronchodilators*
 - a. Go to 'Terms' tab
 - b. Select *Medications* → *Respiratory agents* → *Bronchodilators*
6. Drag and drop
 - a. Drag the *Acute respiratory infections* term into 'Find Patients> and' group panel'
7. Run Query
 - a. Click 'Find Patients> Run Query Options' → 'Run Query'
8. Get Results
 - a. See the 'Query Status' tab for the query result
 - b. The query is now stored in the 'Queries' tab

2. Webclient Layout

The i2b2 Webclient has five main panels, laid out in the following manner.

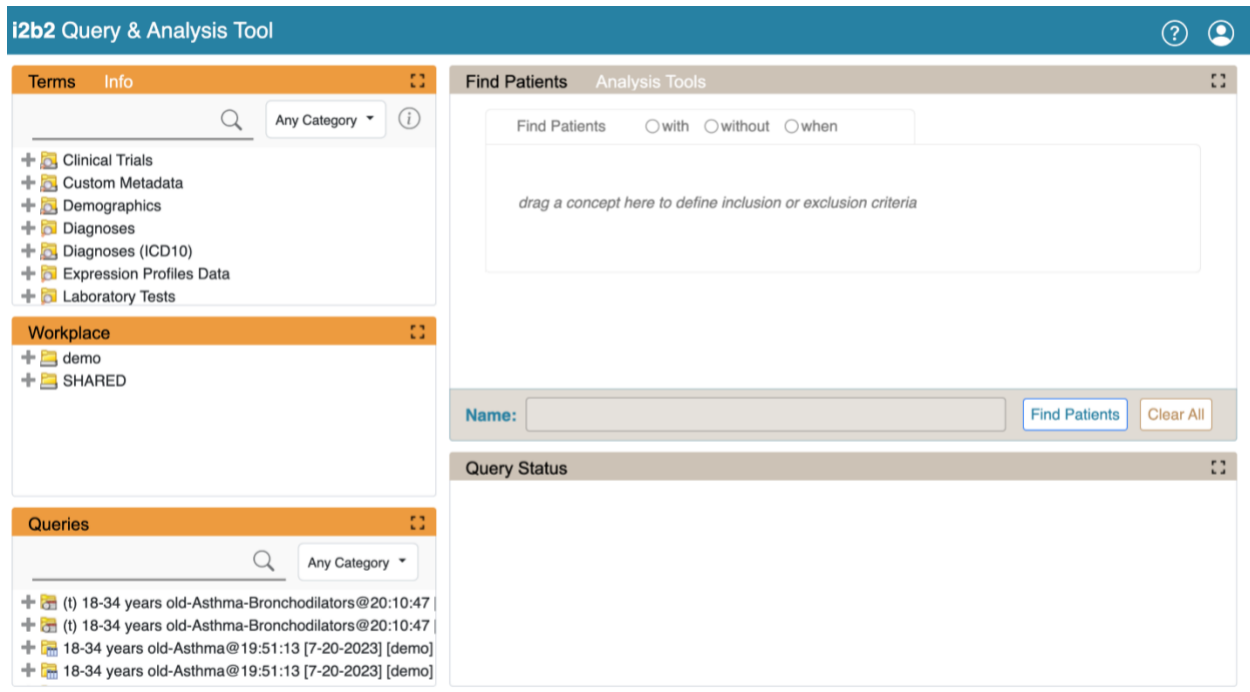

[Find Patients](#) (Top Right)

The **Find Patients** panel is where queries are constructed and run. Users can specify query criteria in the group panels. Query criteria can be any terms contained in the vocabularies (also called ontologies) in the **Terms** section, as well as past queries or query results.

[Terms](#) (Top Left)

The **Terms** panel contains hierarchical vocabularies of terms for constructing queries. Users can navigate through the vocabularies using the default Terms view, or can search by text or codes using the search input field.

[Query Status](#) (Bottom Right)

The **Query Status** panel shows information about the state of a query, such as if it is currently running, or, when it completes, the results. It can also display various graphical representations of the results. The **Query Report** button provides a textual representation of the query

[Queries](#) (Bottom Left)

Previously-run queries and their results are stored in the **Queries** panel. Users may re-use them by dragging them into the Query Name section or groups in the **Find Patients panel** or save them by

dragging them into the **Workplace** panel. The most recently run queries are displayed at the top. Users can use the **search input field** to search for a prior query they created.


Workplace (Middle Left)

The **Workplace** panel stores information, such as prior queries, that can be referred to or used later. Users may store previously-run queries, patient sets, individual patients, and terms in this panel. Users can create folders to organize the workplace. The workplace is private by default, but there is a SHARED folder that everyone has access to. Users can use the SHARED folder to share terms and queries with other colleagues.

3. Query Tool

The **Query Tool** panel is where queries are constructed and run.

3.1 Layout

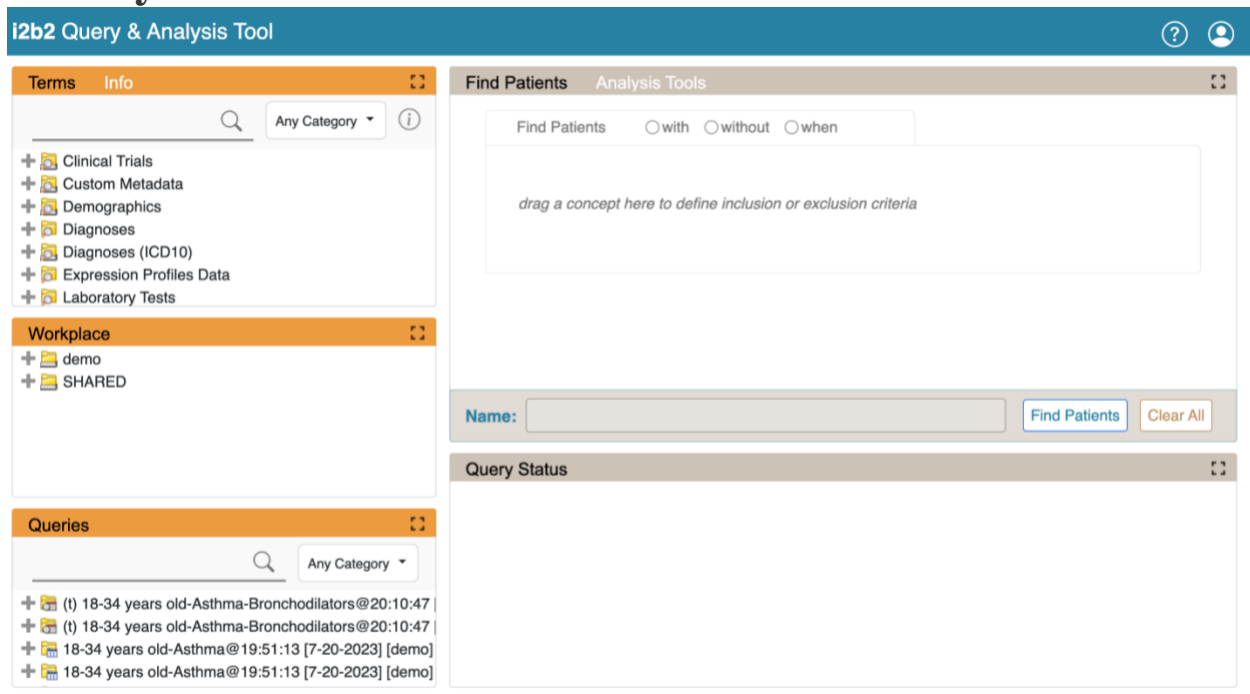

Query Run Options allow users to (1) set options for the panel, and (2) resize the panel.

Query Name shows the name of the current query. Users can also drop a previous query here to load it.

Query Timing Selector enables users to select the appropriate timings for their query, including temporal queries.









Query Groups are where users can drag and drop terms to build queries. There are initially one Group, but more are automatically added.

Tool Bar contains buttons to run queries, clear all the panels, and name the query

3.2 Creating a Query

3.2.1 i2b2 Items for Querying

These items may be dragged into Query Groups to construct queries.

Icons	Names	Description
	Terms	Terms (also called concepts) are the words or phrases that make up vocabularies. Terms are the basic query items.
	Term Groups	Related terms are grouped into hierarchies. When a term group is used in a query, all terms in the group are searched for.
	Modifiers	Modifiers are descriptors that can be attached to terms to make small changes to them. For example, a Diagnosis concept may be additionally described by Modifiers "Principal Diagnosis" or "Secondary diagnosis."
	Modifier Groups	Groups of Modifiers, similar to Term Groups.
	Previous Queries	Previous Queries are queries that have been run before and are stored in the 'Previous Queries' panel.
	Individual Patients	If selected when running a query, a query can return a Patient Set, which is a set of patient numbers that satisfy the query criteria . The Patient Set is stored with the query in the Previous Queries panel. The Patient Set can be expanded, and individual subjects dragged into query groups to be used in a new query.
	Patient Sets	If selected when running a query, a query can return a Patient Set, which is a set of patient numbers that satisfy the query criteria. The Patient Set is stored with the query in the Previous Queries panel and can be dragged into query groups to be used in a new query.
	Encounter Sets	If selected when running a query, a query can return an Encounter Set, which is a set of encounter numbers that satisfy the query criteria. The Encounter Set is stored with the query in the Previous Queries panel and can be dragged into query groups to be used in a new query .

3.2.1 Operators and Order of Precedence

i2b2 items are combined using Boolean operators **AND** and **OR**. Items within each Group are first **OR** ed together. The Groups are then **AND** ed together.

Example 1 (AND): This query finds patients with *Acute Myocardial Infarction AND Angina Pectoris*.

The screenshot shows the i2b2 Query & Analysis Tool interface. On the left, the 'Terms' panel shows a search for 'angina pectoris' with a dropdown menu set to 'Any Category'. Below this, a tree view of medical terms is visible, including 'Diagnoses', 'Circulatory system - 66', 'Ischemic heart disease - 17', and 'Angina pectoris - 3'. The 'Workplace' panel shows 'demo' and 'SHARED' folders. The 'Queries' panel lists several queries, including 'Acute respiratory infections@10:26:30 [7-21-2023] [demo]'. The main 'Find Patients' panel is titled 'Analysis Tools' and contains two query panels. The first panel is for 'Acute myocardial infarction' and the second is for 'Angina pectoris'. Both panels have radio buttons for 'with', 'without', and 'when', with 'with' selected. Below the panels, there are radio buttons for 'and', 'with', 'without', and 'when', with 'and' selected. A 'Name' field contains 'Acute myocardial infarction-Angina pectoris@10:31:18'. A 'Find Patients' button and a 'Clear All' button are visible. The 'Query Status' panel is empty.

Example 2 (OR): This query finds patients with *Acute Myocardial Infarction OR Angina Pectoris*.

The screenshot shows the i2b2 Query & Analysis Tool interface. On the left, the 'Terms' panel shows a search for 'hypertensive disease' with a dropdown menu set to 'Any Category'. Below this, a tree view of medical terms is visible, including 'Hypertensive disease - 40', 'Ischemic heart disease - 17', 'Acute myocardial infarction - 7', and 'Angina pectoris - 3'. The 'Workplace' panel shows 'demo' and 'SHARED' folders. The main 'Find Patients' panel is titled 'Analysis Tools' and contains two query panels. The first panel is for 'Acute myocardial infarction' and the second is for 'Angina pectoris'. Both panels have radio buttons for 'with', 'without', and 'when', with 'with' selected. Below the panels, there are radio buttons for 'and', 'with', 'without', and 'when', with 'and' selected. A 'Name' field contains 'Acute myocardial infarction@8:22:55'. A 'Find Patients' button and a 'Clear All' button are visible. The 'Query Status' panel is empty.

Example 3 (AND of ORs): This query finds patients with (*Acute Myocardial Infarction OR Angina Pectoris*) AND who are male.

*Note that the OR expression in the parentheses is evaluated first.

The screenshot displays the i2b2 Query & Analysis Tool interface. On the left, the 'Terms' panel shows a hierarchical tree of categories: Clinical Trials, Custom Metadata, Demographics (Age, Gender - 133, Income - 134, Language - 133, Marital Status - 133, Race - 134), and Workplace (demo, SHARED). The 'Find Patients' panel on the right shows a query builder with two main panels. The first panel contains 'Acute myocardial infarction' and 'Angina pectoris' connected by an 'or' operator. The second panel contains 'Male' connected by an 'and' operator. The query is named 'Acute myocardial infarction-Male@8:21:33' and includes 'Find Patients' and 'Clear All' buttons.


3.2.2 Constraints

Constraints are used to narrow down the search criteria to a particular time frame, number of occurrences or a specific value or range of values. In i2b2, there are several types of constraints that can be defined in the Find Patients view.


1. Group Constraints

1. Date Constraint

Setting Date Constraint for the entire Group

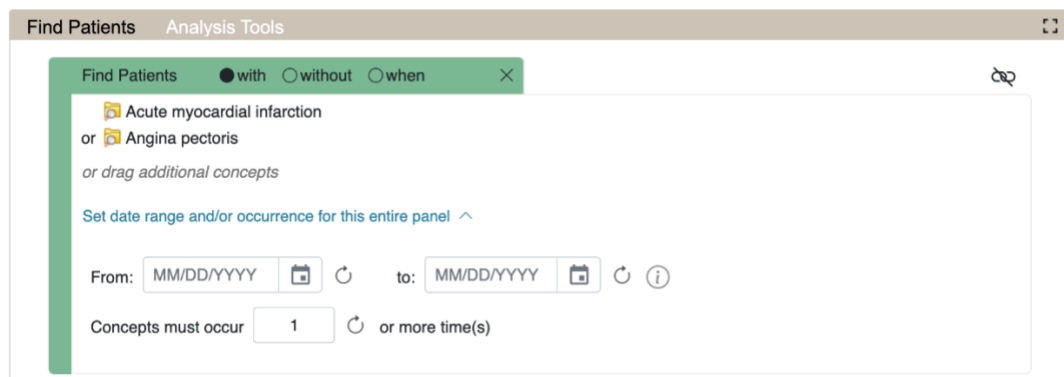

Clicking on 'Set date range and/or occurrence for this entire panel' will expand and display controls to set a date range that can be selected to restrict the occurrences of the items in the Group to within that range. Check one or both of the check boxes to indicate whether there is a begin date (From) or an end date (To) or both.


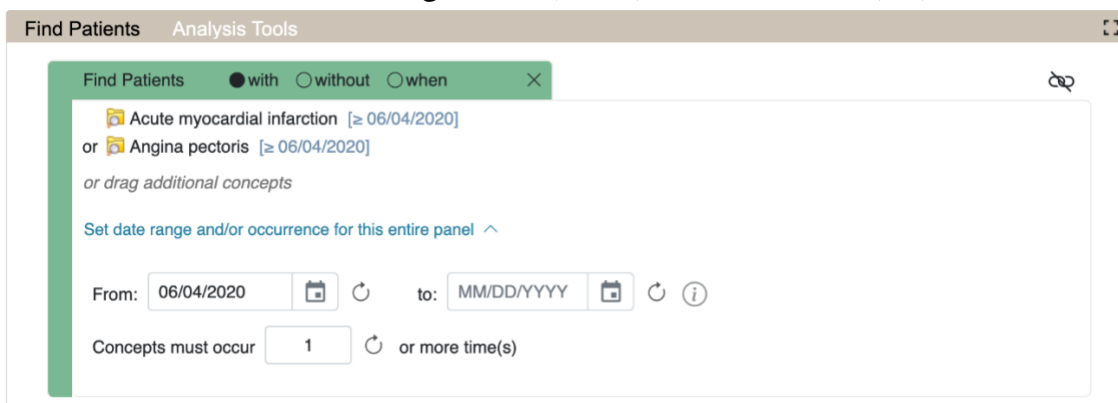


Then type the date directly in the text box or click on the drop-down arrow to open the calendar below to select a date. This will set the date range for all items in this Group.

Setting Date Constraint on individual items

Sometimes only some items in the Group need date constraints or different items in one group require different date constraints. In this case, hover over the item and select **Set Date Constraint calendar icon** to set the date constraint for that item

specifically.

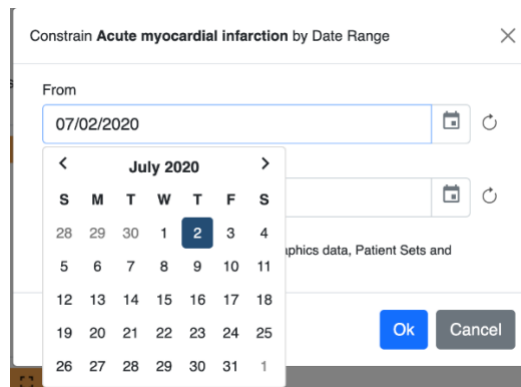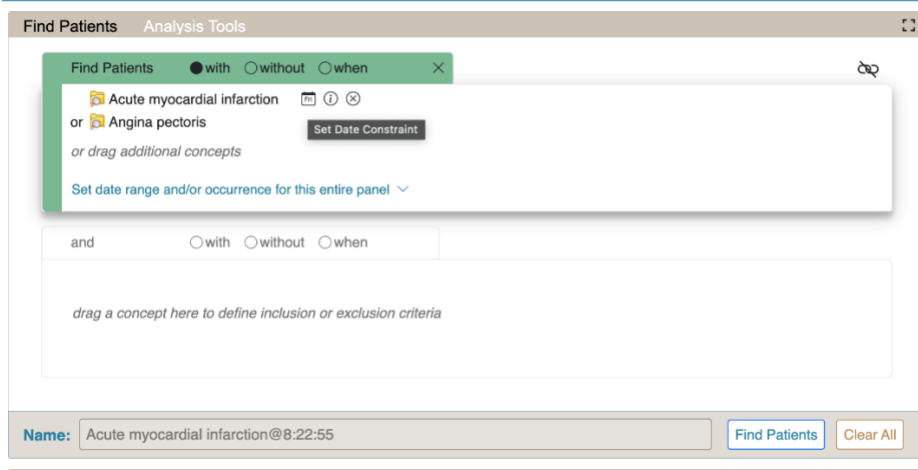


Once a date constraint is applied, **Dates** will now be underlined to indicate a date constraint has been set.




2. Occurrence Constraint

Users can define the number of times an item must occur for each subject.

Clicking on 'Set date range and/or occurrence for this entire panel' will expand and display controls to select to restrict the occurrences of the items in the Group to within that range. Enter a value in the field

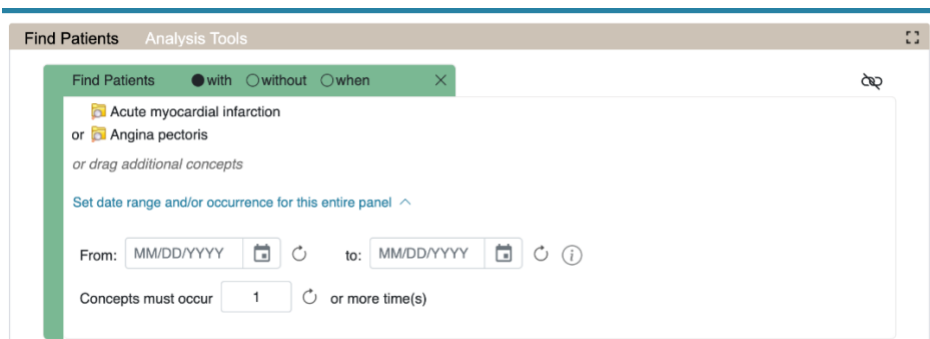


Occurs has been set and the 1 will be replaced with the selected number.

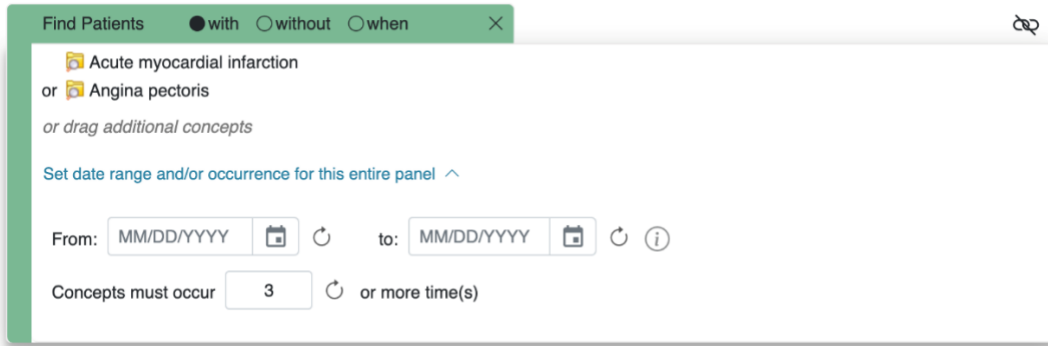

3. Exclusion



Constraint will not include patients with the concepts listed in the panel.

Click on the “without” radio button located in the selection criteria section of the panel. The background color will change to a red hue.

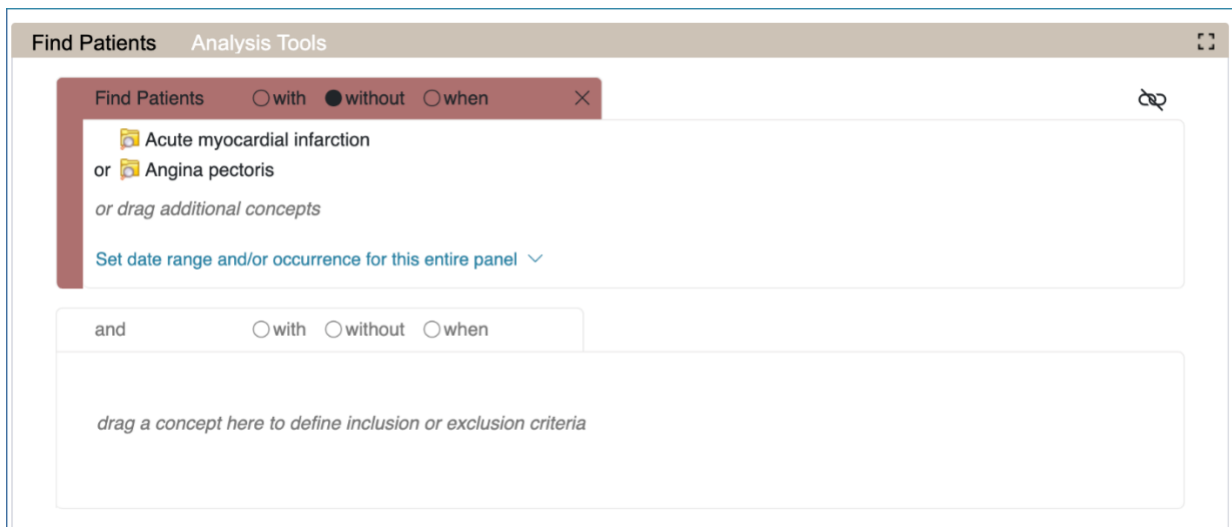

4. Value Constraints (eg. labs)

Constraint is by the value associated to a concept or modifier. Not all items have values. The value constraint window will open automatically when you drop a concept into a panel in the Find Patients view. To set the constraint at a later time, hove over the term and select the corresponding icon to select Set Value.

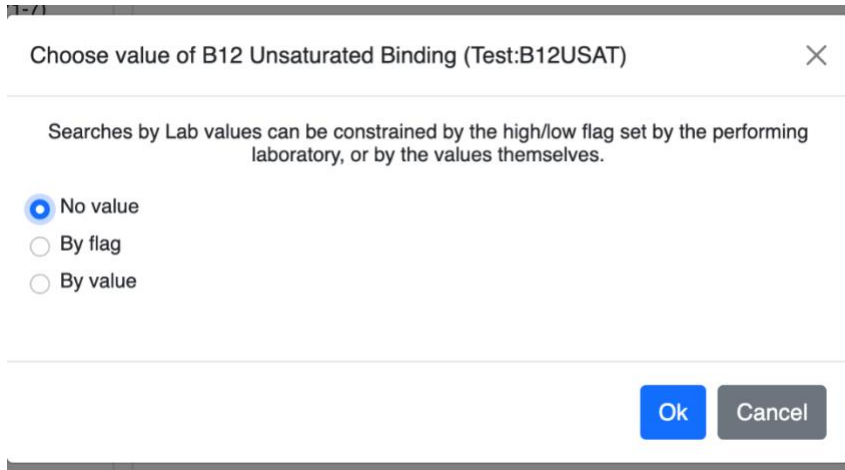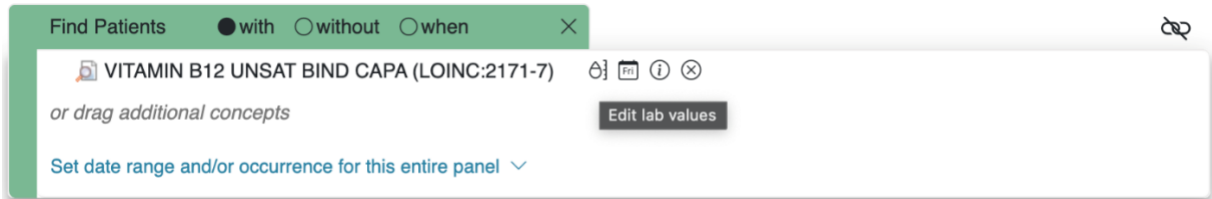

1. Numeric

The numeric value is defined in the source system when the result is entered. Enter a numerical value at the Please enter a value field.

Unlike flag and text values, numeric values use operators, such as greater than or less than when setting the constraints.

The Units will default to the normal units defined for the item. Click on the drop-down arrow to change the units.

(In order for unit conversion to work a parameter called CRC_ENABLE_UNITCD_CONVERSION must be defined for the project. If this is not defined you need to contact your system administrator.)

Choose value of B12 Unsaturated Binding (Test:B12USAT) ✕

Searches by Lab values can be constrained by the high/low flag set by the performing laboratory, or by the values themselves.

No value
 By flag
 By value

Please select operator:
 LESS THAN (<) ▼

Please enter a positive integer value:

Click on a bar segment to help specify a value or range:
 Range in *ng/l*

568 1000 1252 2000

Units = ng/l ▼

2. Categorical

Choose value of Medication Route (Test:MED:ROUTE) ✕

Searches by Lab values can be constrained by the high/low flag set by the performing laboratory, or by the values themselves.

No value
 By value

Please select a value:

inhalation
 injection
 Intravenous
 By Mouth
 By Rectum

3. Text

Choose value of Patient Instructions (Test:MED:INST) ✕

You are allowed to search within the narrative text associated with the term Choose value of Patient Instructions (Test:MED:INST)

No Search Requested Enter Search Text:

Search within Text

Use Database Operators (*Advanced Searching*)

4. Flags (high, low, etc.)

Choose value of B12 Unsaturated Binding (Test:B12USAT) ✕

Searches by Lab values can be constrained by the high/low flag set by the performing laboratory, or by the values themselves.

No value Please select a range:

By flag

By value

5. Modifier Value Constraints

Choose value of Medication Route (Test:MED:ROUTE) ✕

Searches by Lab values can be constrained by the high/low flag set by the performing laboratory, or by the values themselves.

No value
 By value

Please select a value:

- inhalation
- injection
- Intravenous
- By Mouth
- By Rectum

3.2.3 Query Timing


Choose the relationship between the query panels. The default timing is independent.

1. *Independent*

- The items (facts) can occur at any time in the patient's history.

[Example]

Patient has a medication with a dose greater than 500 mg and a medication that is taken once a day documented at any time in their history.

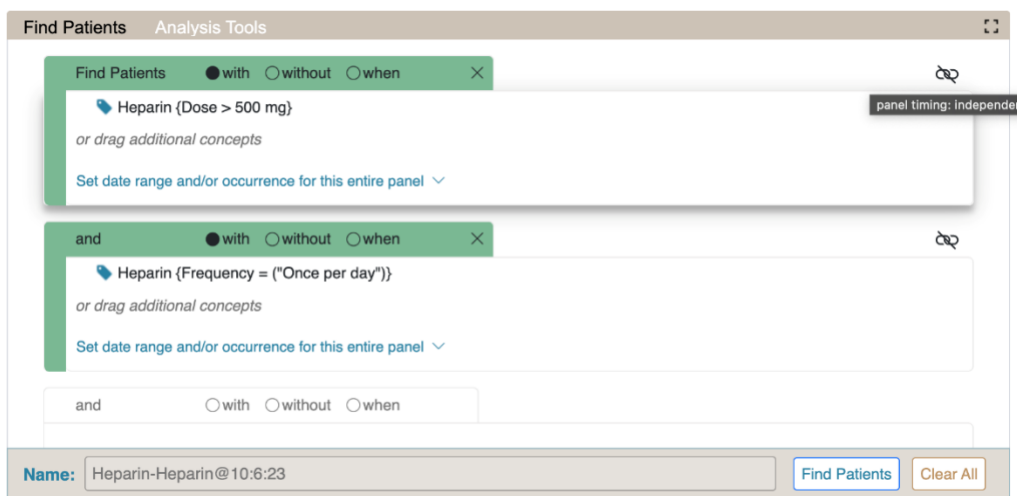



2. Same Financial Encounter

- The items (facts) occur during the same visit (financial encounter).
- When this is selected, individual panels can be included in the "same encounter" query, or they can be treated independently

[Example]


Patient has a diagnosis of Acute MI documented at any time in their history and they have a medication with a dose greater than 500 mg and a medication that is taken once a day documented during the same visit.

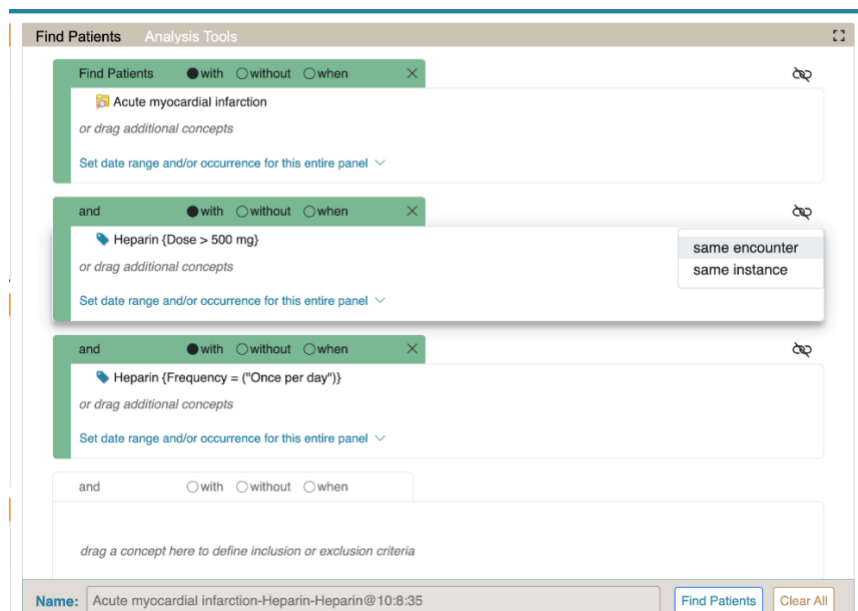

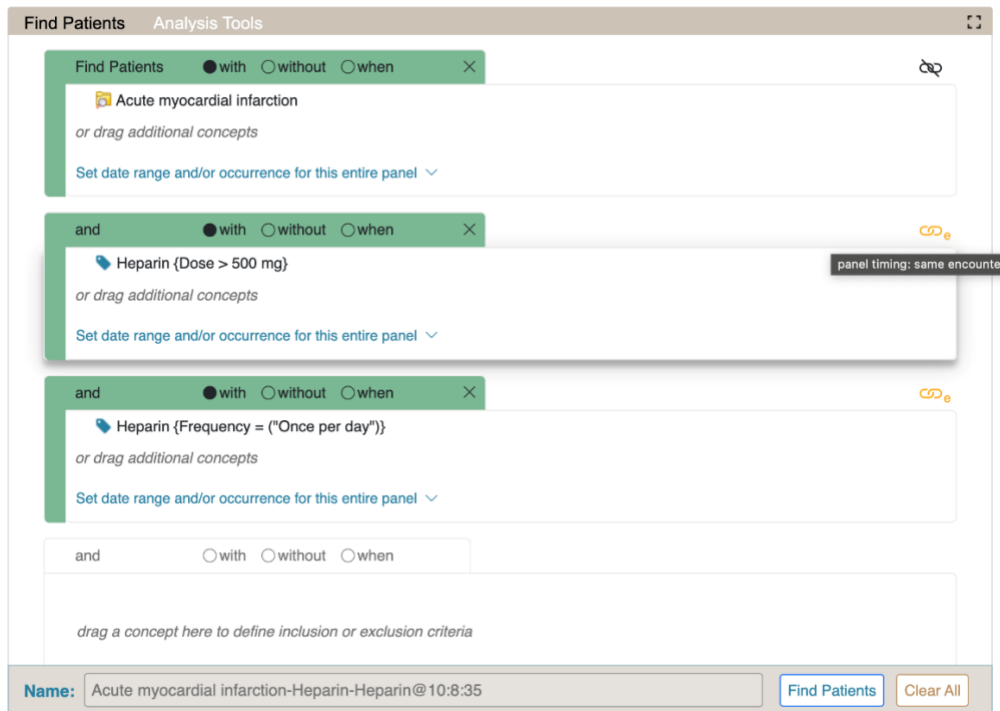


3. Temporal Query

- In this mode, a sequence of events are defined that are anchored by a constraining population. The constraining population is defined exactly like other query types. To build a sequence of events, select the “when” radio button located in the selection criteria section of the panel. The background color will change to a yellow hue.

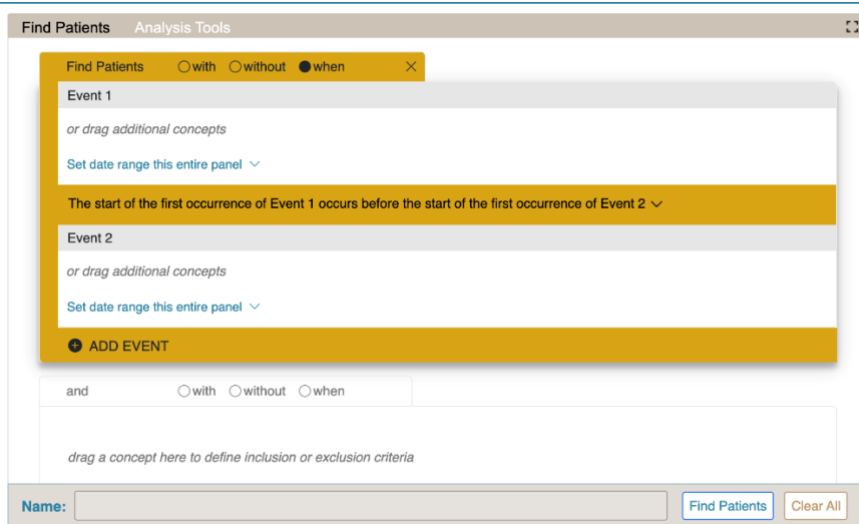

[Example]


- Finds patients who were diagnosed with diabetes (*Event 1*) prior to the first observation of an antidiabetic medication (*Event 2*). The population is constrained to patients over 64 years of age (*bottom panel*).

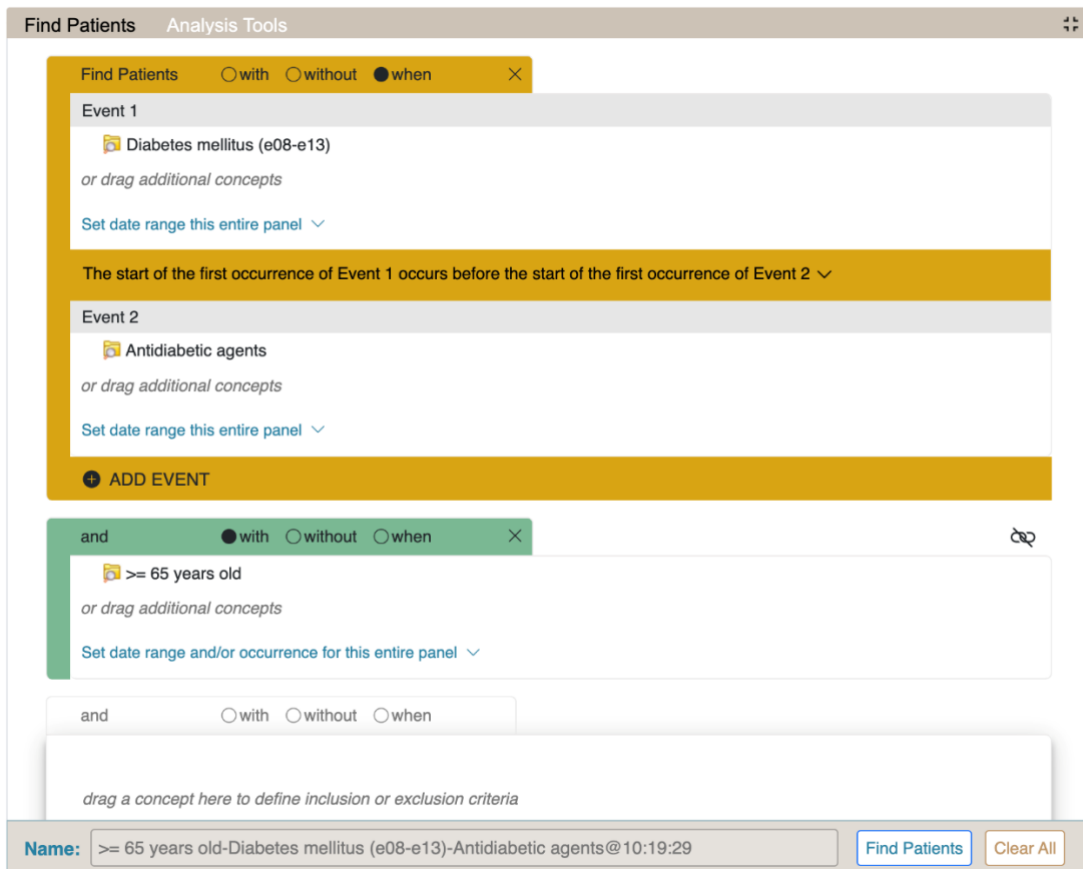


3.2.4 Reloading a Previous Query

Saved queries (either from Queries panel or the Workplace) can be reloaded by dragging them to the Name bar at the bottom of the Find Patients panel. They can be used within a query by dragging them into a query panel. See the [Query section](#).

- Reload a previous query.

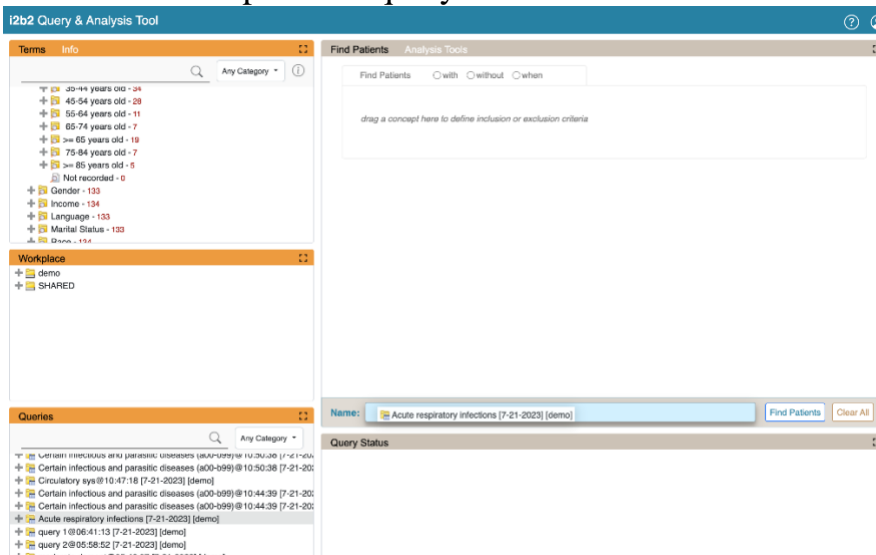

- Use a previous query within a new query.

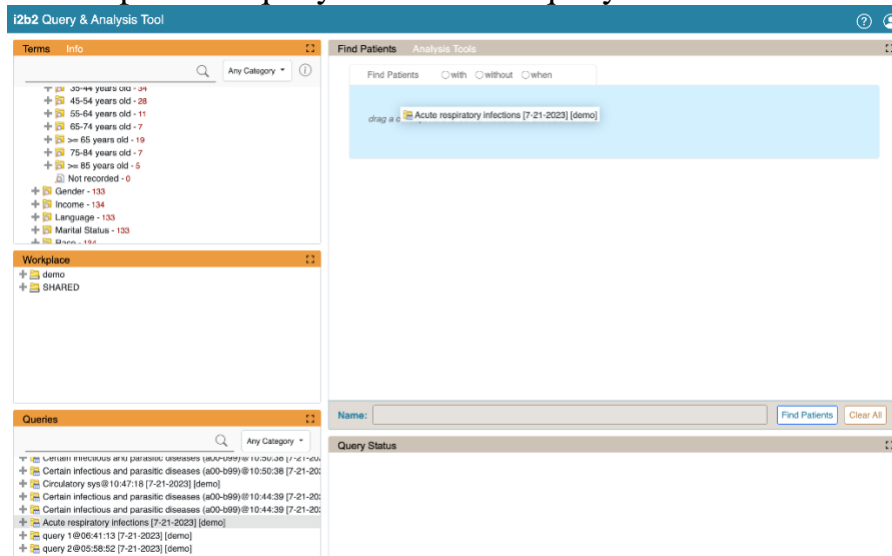


3.2.5 Resetting the Query Tool

There are three options for resetting the **Find Patient panel**. Users can (1) reset all the groups, (2) clear all items in a panel, or (3) delete an individual item from the group.

1. Reset All Groups: Click on the **Clear All button** located on the toolbar at the bottom of the *Find Patients panel*.

2. Clear Items in a Panel / Group : Click on the **Delete “X” button** located in the top right corner of the panel (next to the panel controls).

3. Delete Single Item: Hover over a concept, and click on the X icon to remove from the panel

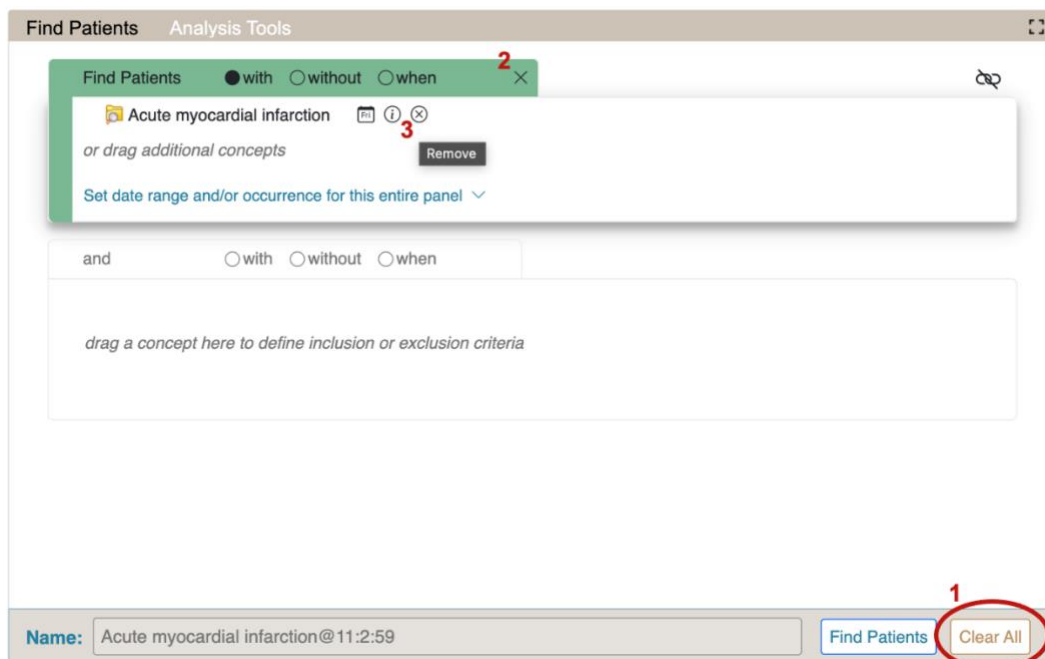

3.3 Running a Query


3.3.1 Run Query

Once the query criterion has been defined, the final steps are to click on the **Find Patient button**, select the result type(s), and click on the 'Run Query' button to run the query.

Once you have defined the search criteria, click on the **‘Find Patients’** button located on the toolbar at the bottom of the Find Patients panel. A new **‘Run Query Options’** window will open.

- A name for the query will default at the **‘Name for query’** field. Either change the name of the query or accept the default.
- In the query result type section, click on the checkboxes for those result types you want to include when the query is run.
- Click on the **‘Run Query’** button.
- The Run Query Options window will close and the querying process will begin.

Once you begin the process of running the query the Query Status view will display "Running Query", the name of the query, and the time that has elapsed.

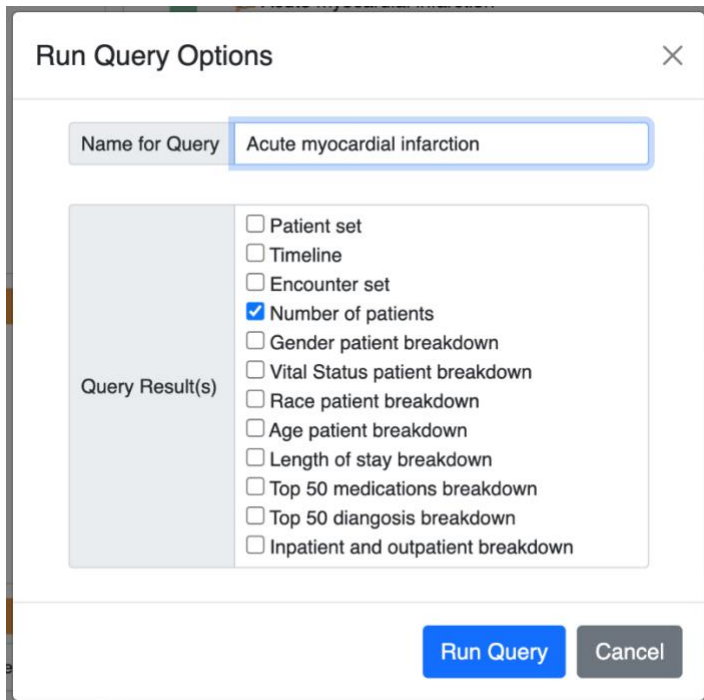

The screenshot shows a dialog box titled "Run Query Options" with a close button (X) in the top right corner. Below the title bar, there is a text input field labeled "Name for Query" containing the text "Acute myocardial infarction". Below this field is a section labeled "Query Result(s)" containing a list of checkboxes. The "Number of patients" checkbox is checked, while all other checkboxes are unchecked. The list of options includes: Patient set, Timeline, Encounter set, Number of patients, Gender patient breakdown, Vital Status patient breakdown, Race patient breakdown, Age patient breakdown, Length of stay breakdown, Top 50 medications breakdown, Top 50 diagnosis breakdown, and Inpatient and outpatient breakdown. At the bottom of the dialog box, there are two buttons: "Run Query" (highlighted in blue) and "Cancel" (greyed out).


3.3.2 Query Result Options

There are three categories of query type: counts, sets, and breakdowns.

Note: The Timeline type is not used in the Web Client. If selected it will work the same as patient set.

1) *Number of Patients*


This option will return a total count for those patients who meet the defined criteria. This is the most common and fastest query type.

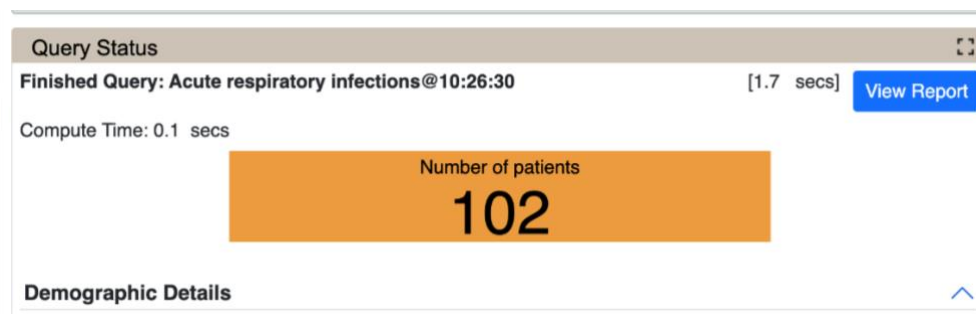

2) *Patient/Encounter Set*

These options return a set of patients or encounters that meet the defined criteria. The set is visible in the Queries panel and can be used in future queries, the workspace, and by plugins.

Only information on the query's completion is shown in the Query Status panel.

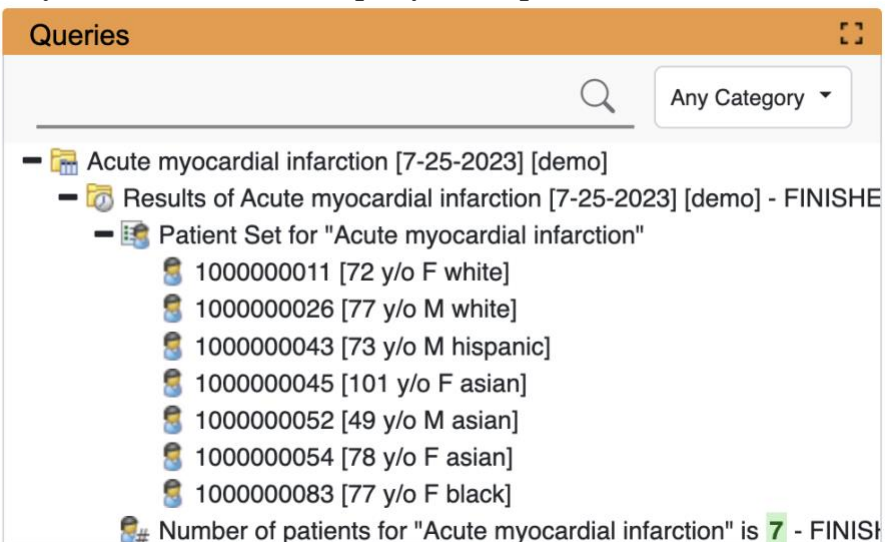

3) Patient Breakdowns

The following query types will return patient counts stratified by the selected breakdown category (e.g., Gender or Age).

Breakdown type

Gender patient breakdown

Vital Status patient breakdown

Race patient breakdown

Age patient breakdown

Length of Stay breakdown

Top 20 Medications Breakdown

Top 20 diagnoses Breakdown

Inpatient and Outpatient Breakdown

3.3.3. Query Status

When a query is running, the query status panel will show its status.

Status	Meaning
---------------	----------------

Running	The query is currently running on the server and the client is waiting for results.
---------	---

Finished	The query has completed. Results will be shown below this message.
----------	--

Queued	The query took longer than the configured timeout. (Defaults to 180 seconds.) It is still running in the background and the results will appear in Previous Queries when it completes.
--------	--

Error An error occurred running the query. More information might be available in the XML Message Log.

4. Navigate Terms

Terms, (also called concepts), make up vocabularies (or ontologies) in i2b2 and are the building blocks of i2b2 queries. The Terms Panel is used to find terms for constructing queries. It has two modes, **Term view** and the **Search** view.


The **Terms** tab is used to drill into a hierarchy of terms and find terms of interest to query.


The **Search input** is used to search for terms or codes to be used in queries.


Navigate Terms Tab

The **Navigate Terms** tab is used to find ontology items to be used in a query. This panel displays the i2b2 query terms or concepts. Terms are grouped with similar items and displayed in a hierarchical tree structure that can be expanded to view all items in the category.

Folders are expanded by clicking on the '+' sign. The '+' changes to a '-' and the next level of detail is displayed. The hierarchical structure consists of [containers](#), [folders](#) and [leaves](#).

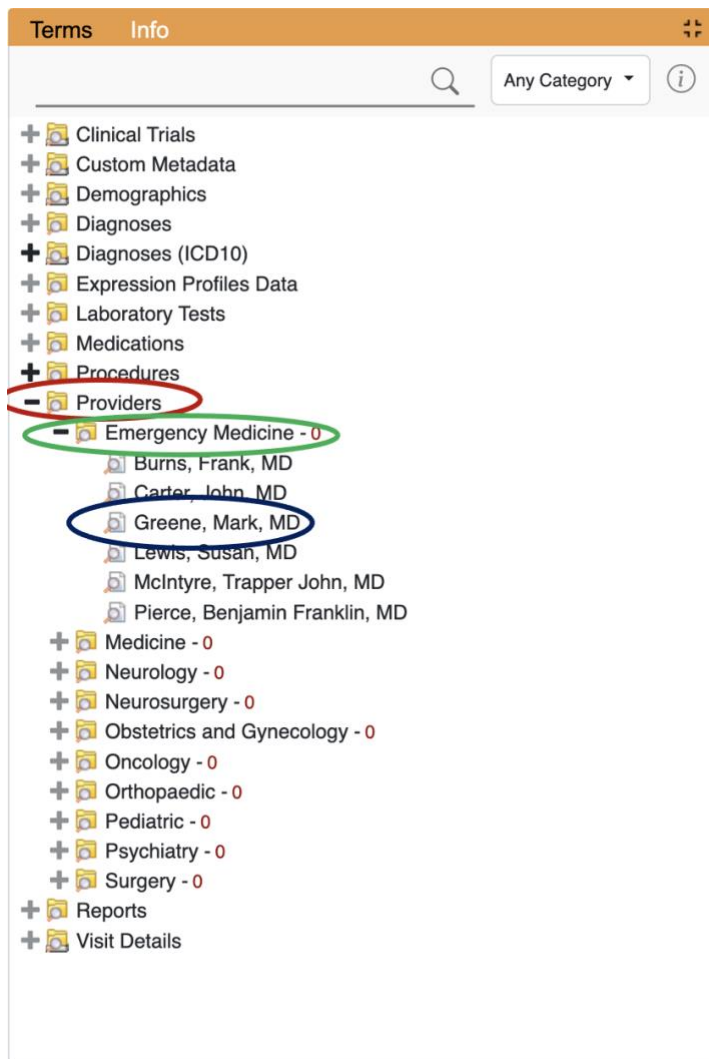

 **Containers** are the highest level of the hierarchy. They can be expanded by clicking on the *plus sign* (+) or collapsed by clicking on the *minus sign* (-). Containers cannot be dragged to other panels.


 **Folders** Like containers, folders can be expanded using the + and collapsed using the -. Unlike containers, folders can be dragged to other panels. When a folder is dragged into a query group, the resulting query will return results for everything in the folder.



Leaves are the lowest level of the hierarchy. They cannot be expanded or collapsed. They can be dragged into query groups as query criteria.

Find Mode


The **Find** mode provides a way to search for a query items. Query terms can be located using the Find input field.

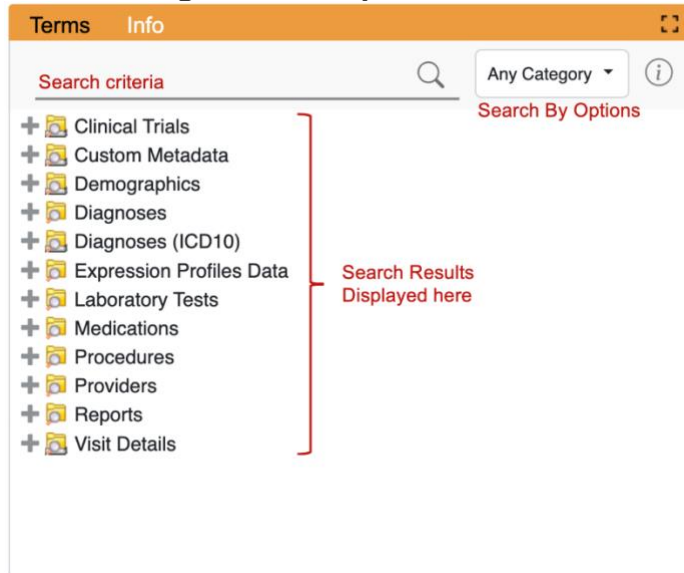

A. The Search by Names mode is used to find all query items containing particular text.

1. Enter search string. All results returned will be those terms whose name or code CONTAINS the text entered in the search input field.
2. Choose a category from the drop-down menu or leave it as the default of 'Any Category', which will search all categories.
3. Click Enter on the keyboard or click the search icon to run the search
4. Select any terms returned by the search and drag them into a query group to use as query criteria.



[Example] Find all terms containing the text 'rheumatic fever'.

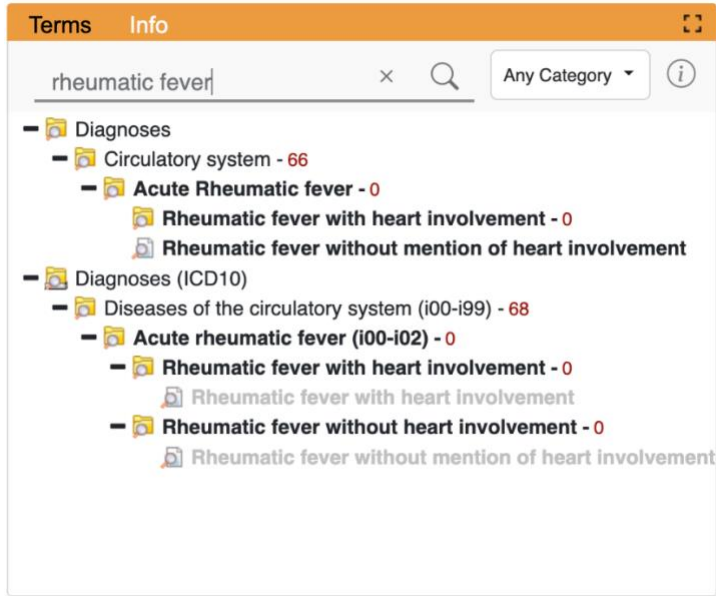

Note: New in v1.7.12

The Find Terms view has changed. The same query now shows many fewer results and it shows their placement in the hierarchy. Search results are highlighted in bold emphasis, and parent results can also be dragged into other panels. Results which are hierarchically within another result are omitted to save space (e.g., many of the results in the graphic above are within the Acute Rheumatic fever folder).

Right click on a term to view a menu. Click “View in Tree” to restore the navigate mode.

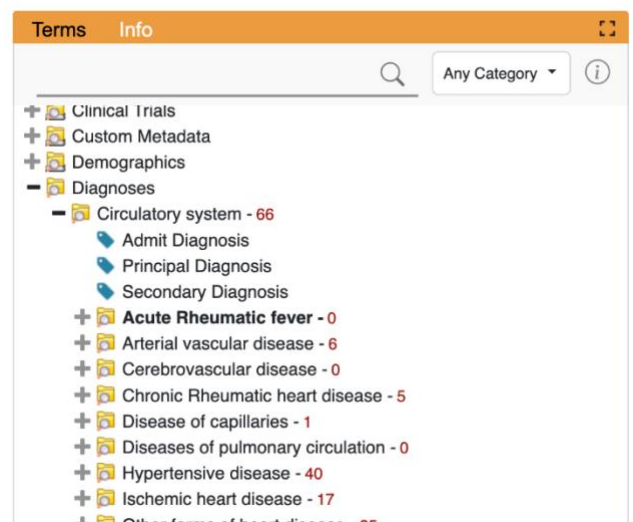
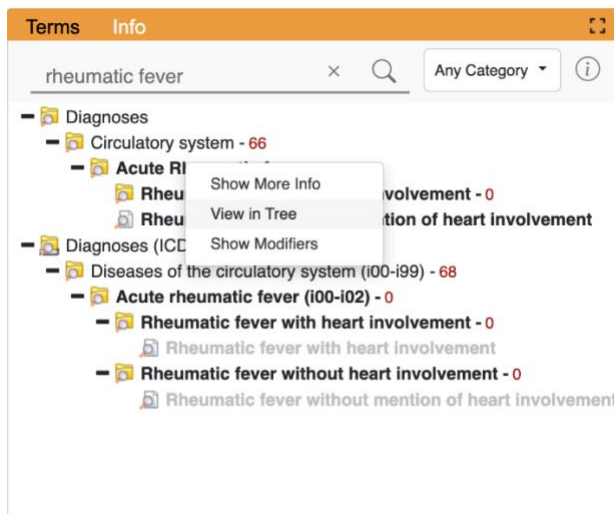


B. The Search by Codes tab is used to find coded variables.

To use this option, a code must be entered and a coding system selected.


[Example] Find LOINC code 32673-6

1. Enter the code in the search input field
2. Click on the **drop-down arrow**, select **Coding System**, and select **LOINC** from the drop down list.
3. Click Enter on the keyboard or click the search icon to run the search
4. The results will be displayed in the results section.

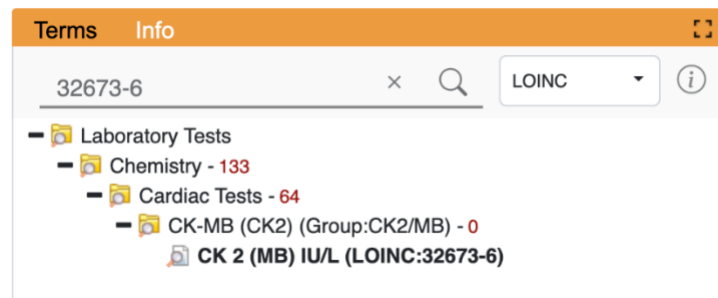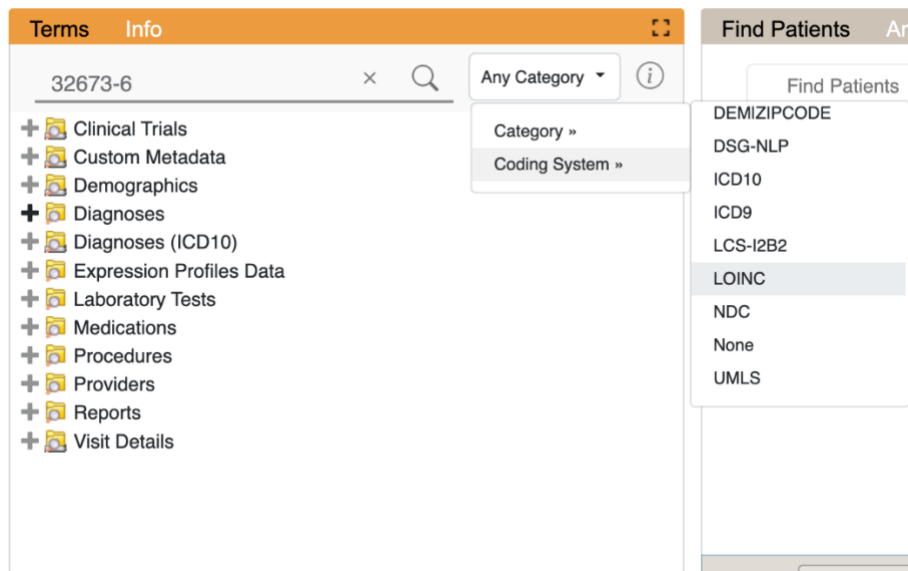

Options for Finding Terms

Right clicking on Terms opens an the Options for Navigating/Searching Terms that provides additional options for finding terms.

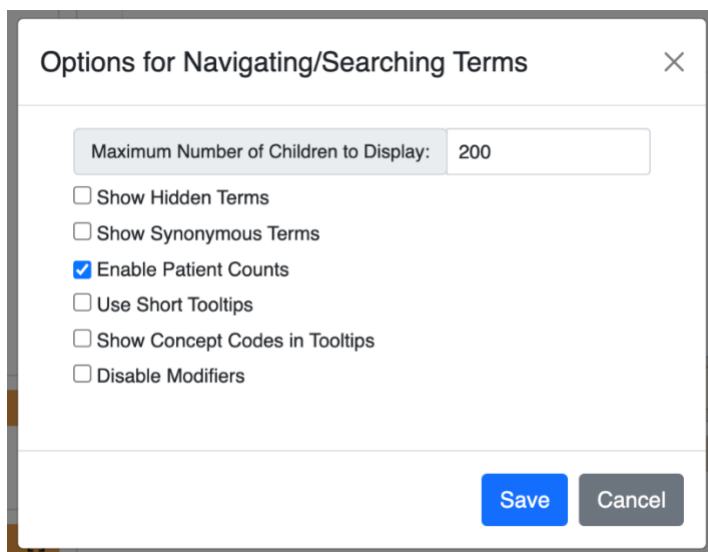

Maximum Number of Children to Display

The default value for **Maximum number of children to display** for a search of terms or codes is usually 200. However, this number can be changed to retrieve either more or fewer results. When a search is run, once the maximum number of terms is reached a message will appear to say that the query has reached the maximum number of terms and suggest searching again with a more specific query. Running a search with a more specific query will narrow results as the selection criteria will not be as broadly defined.

Show Hidden Terms

In the metadata table, an attribute is set to determine whether a term is active, inactive or hidden. If an item is marked as hidden it will not appear in the list of terms. Selecting **Show Hidden Terms** will display those items marked to be hidden. For easy identification, hidden terms will be **red** when displayed in the results list.

Show Synonymous Terms

Some terms may be flagged in the metadata table as a synonym of another term. Selecting **Show Synonymous Terms** will display both the original term and the synonym in the *Terms view*. For easy identification, synonymous terms will be **blue** when displayed in the results list.

Enable Patient Counts

The **Enable patient counts** option will display the total number of patients that have a fact entered for that item.

Use Short Tooltips

In the i2b2 webclient, a tooltip displays when the cursor hovers over an item. By default, the tooltip that displays contains the entire path name and the tooltip that is specific to that item. Enabling the **Use short tooltips** option will only display the tooltip information; the path will not be included.

Show Concept Codes in Tooltips

In addition to showing the tooltip, this option will also display the concept code associated with the item.

Disable Modifiers

The ability to search for modifiers will not be available if the **Disable modifiers** option is checked off.

5. Query Status


The **Query Status** panel contains information regarding the state of a query.

Show Query Status

While a query is running, the name of the query and how long it has been running will appear.

Additional information that shows up in this tab depends on which options are selected in the 'Run Query Options' popup box when the query is run.

Once the query has finished, the number of patients meeting the search criteria, and the time it took to run the query, will appear.

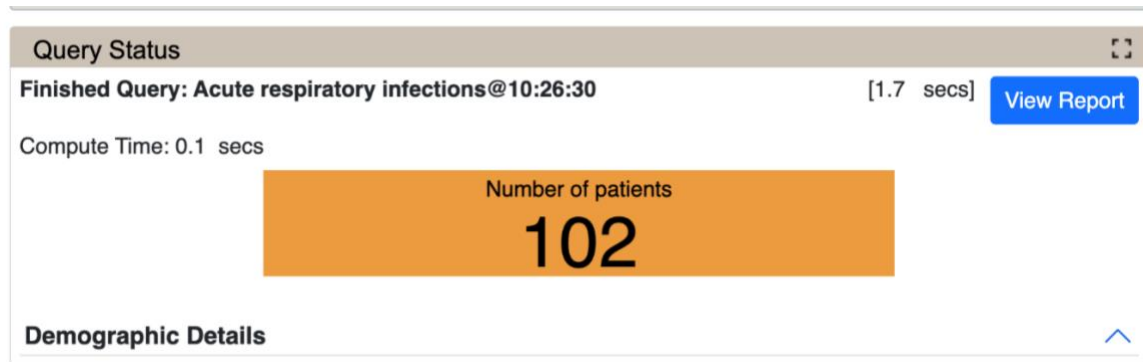

Graph Results

If any of the breakdown categories (Age, Race, Gender, Vital Status, Length of Stay, Top 20 diagnoses, Top 20 medications, Inpatient and Outpatient) are selected in 'Run Query Options', a bar graph will display for each category selected.

This is the Run Query Options showing default breakdown types:


Run Query Options

Name for Query: Acute respiratory infections

Query Result(s)


- Patient set
- Timeline
- Encounter set
- Number of patients
- Gender patient breakdown
- Vital Status patient breakdown
- Race patient breakdown
- Age patient breakdown
- Length of stay breakdown
- Top 50 medications breakdown
- Top 50 diagnosis breakdown
- Inpatient and outpatient breakdown

[Run Query](#) [Cancel](#)

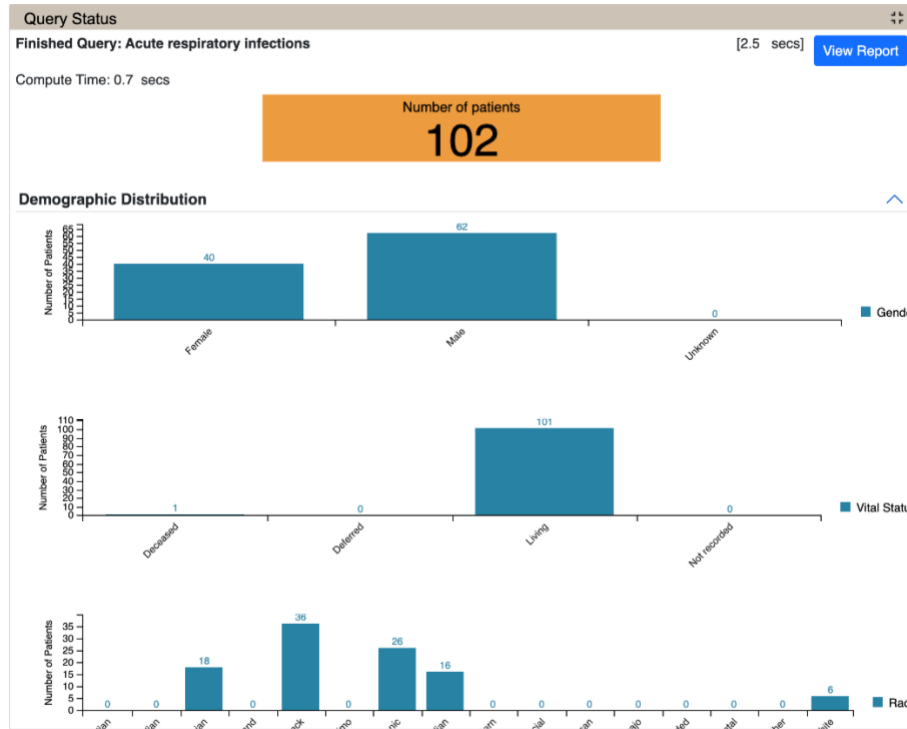

Query Report

The 'Query Report' button displays a textual representation of the query. It can be printed by clicking on the 'View Report' button in the upper right hand corner.

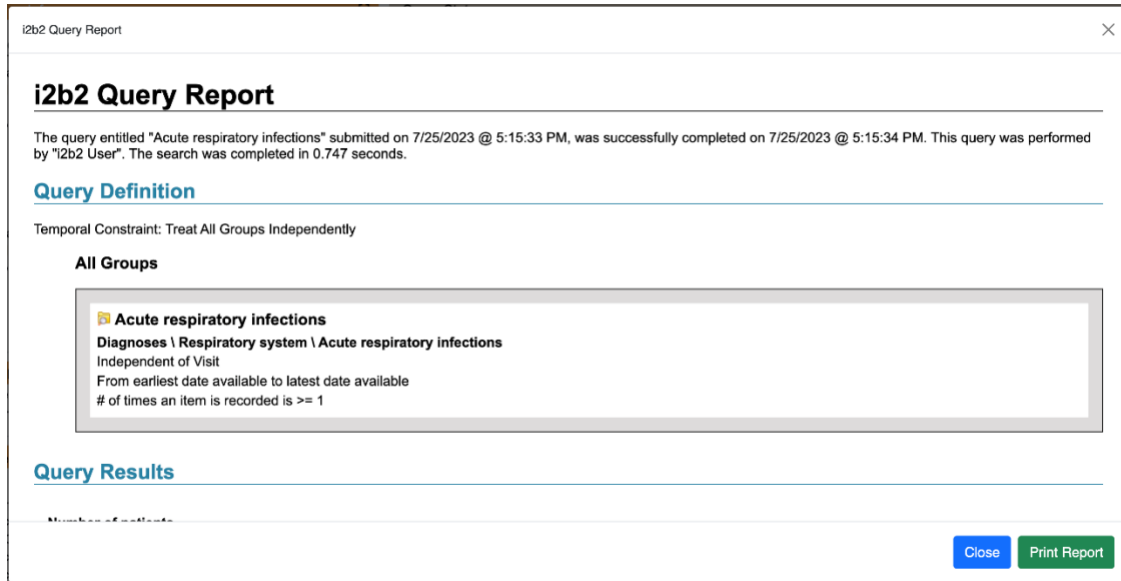

The screenshot shows a window titled "i2b2 Query Report" with a close button in the top right corner. The main content area is titled "i2b2 Query Report" and contains the following text:

The query entitled "Acute respiratory infections" submitted on 7/25/2023 @ 5:15:33 PM, was successfully completed on 7/25/2023 @ 5:15:34 PM. This query was performed by "i2b2 User". The search was completed in 0.747 seconds.

Query Definition

Temporal Constraint: Treat All Groups Independently

All Groups

Acute respiratory infections
Diagnoses \ Respiratory system \ Acute respiratory infections
Independent of Visit
From earliest date available to latest date available
of times an item is recorded is >= 1

Query Results

Number of Patients



At the bottom right of the window, there are two buttons: "Close" (blue) and "Print Report" (green).


6. Previous Queries

View and search queries run by the user.

The 'Queries' panel has 2 views. The default mode lists all queries run by the user, and the search mode allows users to search for their previously run queries.

6.1 'Queries' tab

The 'Previous Queries' tab displays queries run by a user, with the most recent ones listed first. Results associated with the query can also be viewed in this tab. Both a *standard i2b2 query* () and a *Temporal query* () will appear in the 'Queries.' Please see [3. 'Query Tool'](#) for query types.

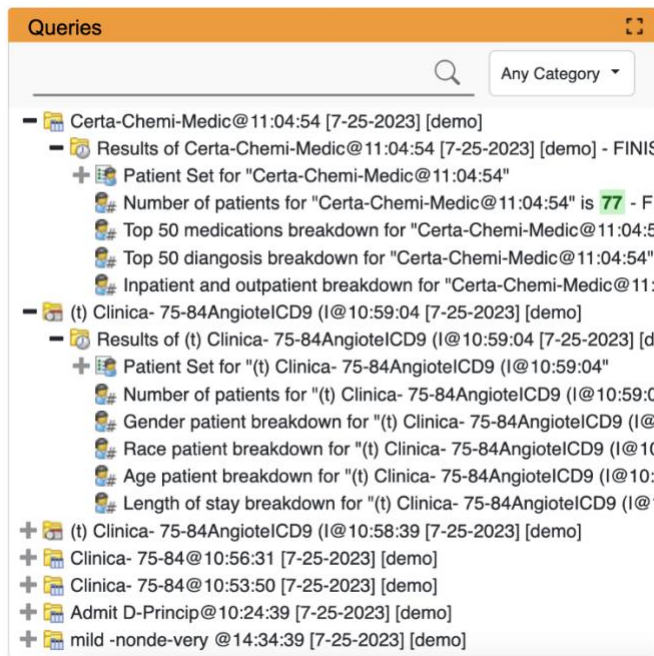


6.2 Query Search mode



The **Find** mode provides a way to search for a query or query result type. Query terms can be located using the Find input field and will find all query items containing particular text.

5. Enter search string. All results returned will be those terms whose name or code CONTAINS the text entered in the search input field.

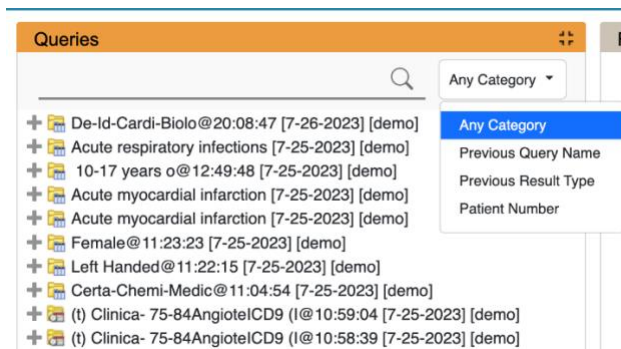6. Choose a result type from the drop-down menu or leave it as the default of 'Any Category', which will search all categories. Other categories are: 'Previous Query Name', 'Previous Result Type', 'Patient Number'.
7. Click Enter on the keyboard or click the search icon to run the search
8. Select any terms returned by the search and drag them into a query group to use as query criteria.



6.3 Layout of a Previous Query Display


A **previous query** ( / ) is created when a query is run from the 'Query Tool'. These previous queries can be used to run new queries or to easily populate other panels.

The query in the 'Previous Queries' panel is comprised of multiple hierarchical levels that can be seen by expanding each level of the tree. These levels are outlined below.



6.3.1 Previous Query Name

The name of the previous query is the top level of the hierarchal tree and is comprised of these three components.

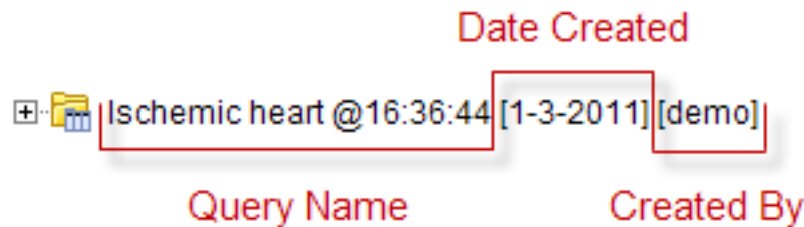1. **Query Name** - the name given at the time the query was created. The default name contains an abbreviated description and the time the query was run.
2. **Date** - the date the query was created.
3. **User Id** - the user who created (ran) the query.




6.3.2 Results Folder

The **Results folder** () contains the results of the query. To see the query results simply expand the folder by clicking on the *plus sign* () located next to the result name.


6.3.2.1 Results

The results associated with a **previous query** are determined by the **query result type(s)** that are selected when the query is run in the 'Query Tool'.


6.3.2.1.1 Patient Set

The **Patient set** () appears if *Patient set* or *Timeline* was checked at the time the query was run. The set contains all the patients who met the defined search criteria.


6.3.2.1.2 Encounter Set

The **Encounter Set** () appears if *Encounter set* was checked at the time the query was run. The set contains all the patients whose encounters met the defined search criteria.

6.3.2.1.3 Number of Patients

The **Number of Patients** () appears if *Number of Patients* was checked at the time the query was run. The results are a total count for those patients who met the defined criteria.


6.3.2.1.4 Breakdown

The **breakdown** () appears if *Gender / Vital Status / Race / Age patient breakdown* was checked at the time the query was run. The results are a breakdown by gender / vital status / race / age for those patients who met the defined criteria.

6.4 Navigation

6.4.1 Menu

Using the **mouse right click** on the 'Queries' panel, users can open the **Queries Menu**. 'Display,' 'Rename,' 'Delete' options act on the selected, or mouse hovered previous query. When the menu is called from a blank space, those menu options are applied to the previous target query.

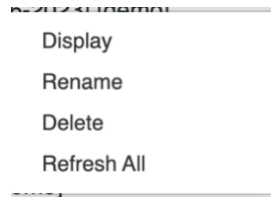

6.4.1.1 Display

Users can see a previous query information by selecting **Display** from the *popup menu*. It loads its query name, the query design, and the results in the 'Find Patients' panel.

6.4.1.2 Rename a Previous Query

Users can rename a previous query from the list by selecting **Rename** from the *popup menu*.

Note: When dragging a previous query to another view the new name will appear when it is dropped

6.4.1.3 Delete / Remove a Previous Query

Users can rename a previous query from the list by selecting **Delete** from the *pop-up menu*.


WARNING: *Once you click on Delete, the previous query will be removed from the view and there is no way to restore it. Therefore, you need to make sure you want to remove it from the list.*

6.4.1.4 Refresh the List

Users can refresh the list of previous queries from the list by selecting **Refresh All** from the *pop-up menu*.

6.4.2 Options

Right click on the 'Queries' label to open an **Options** window. Users can define how many previous queries to display, the order in which they will appear in the list, and auto refresh setting.

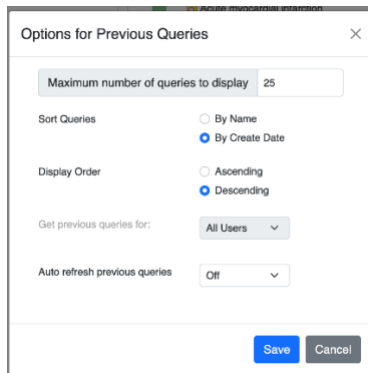

6.4.2.1 Maximum to Display

The **Maximum number of queries to display** defines how many previous queries will be visible in the 'Queries' panel. Once the maximum number is reached the oldest query will no longer appear in the view.

TIP: Recommendation: Use the 'Workplace' panel to save those previous queries that will be used often. In addition to being more convenient it will also prevent you from losing a common query if the maximum number to display has been reached.

6.4.2.2 Sort Queries by Name

Selecting to sort **By Query Name** will sort the queries in the 'Queries' alphabetically by the name of the query. Once you select to sort by name, you need to define whether it will be in ascending or descending order.

Ascending: queries are displayed in alphabetical order with the *beginning* of the *alphabet* at the *beginning* of the list. (A to Z)

Descending: queries are displayed in alphabetical order with the *beginning* of the *alphabet* at the *end* of the list. (Z to A)

6.4.2.3 Sort by Create Date

Selecting to sort **By Create Date** will sort the queries in the 'Queries' by the date and time the query was created. Once you select to sort by the creation date, you need to define whether it will be in ascending or descending order.


Ascending: displays the *oldest* query *first*.

Descending: displays the *newest* query *first*.

NOTE: Renaming a query will not affect the sort order as the date and time is stored with the previous query.

6.4.2.4 Auto Refresh

Users can turn off the auto refresh, or select the auto refresh frequency from the combo box.



NOTE: If the auto refresh option is turned off, users need to use the refresh button () to see the updated previous queries list.

6.5 Using a Previous Query

6.5.1 Draggable Items


Different parts of the previous query can be dragged and dropped (copied) from the **'Queries' panel** to other views in the i2b2 Web Client.

6.5.1.1 Previous Query

Associated with each **previous query** ( / ) are the items used to run the query, any constraints that were defined, and the results. A previous query can be added to other views by dragging the previous query name from **'Queries' panel** to one of the following views.


- **'Find Patients' panel** by drop in the *panel (as a concept)* or *Name section (to reload the query)*.
- **'Workplace' panel** by drop in a *user* or *shared folder*.

6.5.1.2 Patient Set

A specific set of patients can be added to other views by dragging the **Patient set** () from **'Queries'** to one of the following views.

- **'Find Patients' panel** by drop in the panel.
- **'Workplace' panel** by drop in a *user* or *shared folder*.

6.5.1.3 Encounter Set

A specific set of encounters (visit list) can be added to other views by dragging the **Encounter Set** () from **'Queries'** to one of the following views.

- **'Find Patients' panel** by drop in the panel.
- **'Workplace' panel** by drop in a *user* or *shared folder*.

6.5.2 Add to 'Find Patients'

The 'Find Patients' is designed to simplify the process of retrieving information from the database associated to the i2b2 Web Client.



Within the '**Queries**' panel are several items that can be used in the 'Find Patients' to run a query. The following sections explain how to add the previous query, patient set and encounter set to the 'Find Patients'.

6.5.2.1 Add Previous Query to 'Find Patients'

Within the 'Find Patients', a previous query can be added to one of the following two locations.



Name field	The items and constraints from the original query will be used to create a new query. This information can be edited to make a new query or can be used as is to run the query again.
Panel (Group)	Used within a new query (query-in-query). The information associated with the previous query can not be edited. Add additional search criteria (items) to be used in the new query.

6.5.2.1.1 Add Previous Query to the Name field (New Query)




1. Highlight the **previous query** ( / ) by clicking on the name of the query. While holding the *left mouse button* down, drag the item to the 'Query Tool'.
2. Drop the item into the **Name field**.
3. The item(s) associated with the previous query will display in the appropriate group(s), the name of the query will appear at **Query Name**, and any constraints that were defined when the original query was run will now default with the new query.

NOTE: When you run the query it will receive a new name. This is to distinguish the first previous query from the second.




6.5.2.1.2 Add Previous Query to the Panel (Query-in-Query)

1. Highlight the **previous query** ( / ) by clicking on the name of the query. While holding the *left mouse button* down, drag the item to the 'Query Tool'.
2. Drop the item into the **Find Patients panel**.

6.5.2.2 Add Patient Set to the 'Query Tool'

1. Click on the plus sign () next to the name of the previous query that contains the patient set.
2. Click on the plus sign () next to the results folder.
3. Highlight the **Patient set** () by clicking on its name.
While holding the *left mouse button* down, drag the item over to the 'Query Tool'
4. Drop the item into the **Find Patients panel**.

6.5.2.3 Add Encounter Set to the 'Query Tool'



1. Click on the plus sign () next to the name of the previous query that contains the encounter set.
2. Click on the plus sign () next to the results folder.
3. Highlight the **Encounter set** () by clicking on its name.
While holding the *left mouse button* down, drag the item over to the 'Query Tool'
4. Drop the item into the **Find Patients panel**.

6.5.3 Add to 'Workplace' panel




Information in the workplace is related to the most common concepts and queries that an individual uses and in essence becomes their personal workplace.



6.5.3.1 Add Previous Query to 'Workplace' Panel




1. Highlight the **previous query** ( / ) by clicking on the name of the query.
2. While holding the *left mouse button* down, drag the previous query name to the **'Workplace' panel**.
3. Drop the previous query into either your workplace folder or the shared folder.
4. The previous query will now appear in the workplace with the same name.

6.5.3.2 Add Patient Set to 'Workplace' panel

1. Click on the plus sign () next to the name of the previous query that contains the patient set.
2. Click on the plus sign () next to the results folder.
3. Highlight the **Patient set** () by clicking on its name.

4. While holding the *left mouse button* down, drag the **Patient set** to the '**Workplace**' **panel** .
5. Drop the Patient set into either your workplace folder or the shared folder.
6. The Patient set will now appear in the workplace.

6.5.3.3 Add Encounter Set to 'Workplace' panel

1. Click on the plus sign () next to the name of the previous query that contains the encounter set.
2. Click on the plus sign () next to the results folder.
3. Highlight the **Encounter Set** () by clicking on its name.
4. While holding the *left mouse button* down, drag the **Encounter set** to the '**Workplace**' **panel**.
5. Drop the encounter set into either your workplace folder or the shared folder.
6. The encounter set will now appear in the workplace.

7. Workplace


A users' personal workspace, for storing favorite concepts and queries.
The main objectives of the 'Workplace' are:

- Save and organize user specific items.
- Share information with other project users.
- Project managers can view the workplace for individual team members.


The items stored in the 'Workplace' panel can be used in many of the i2b2 panels and plug-ins, such as the 'Find Patients' panel.

7.1 Layout of the 'Workplace' panel

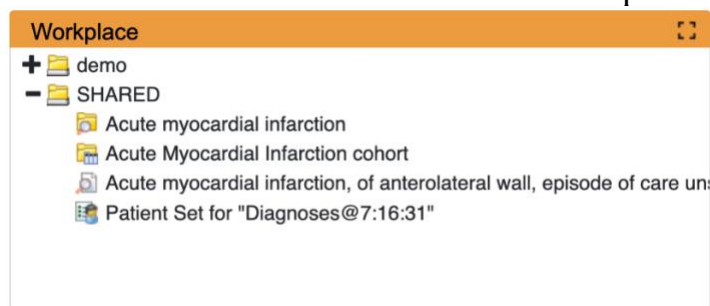


The 'Workplace' panel is hierarchical display of items (concepts, queries) added by the user. There are three main branches or components of the tree that displays.



7.1.1 Containers








Containers are the highest level folders () that are created automatically. The **User container** is created with the user's name and is used to store items added by the user. The **SHARED container** is a place to share concepts and queries with other users. Containers cannot be dragged into 'Find Patients' groups as query criteria.

7.1.2 Folders

A **folder** () contains other folders and / or items added by the user. It can be expanded and collapsed by clicking on the *plus sign* () or *negative sign* () next to it.

7.1.3 Leaves (Items)





Leaves (items) are the lowest level of the hierarchy. They cannot be expanded further. Some examples of items are listed below.

- Concepts ()
- Group of concepts (ontology folder) ()
- Modifiers ()
- Group of modifiers (modifier folder) ()
- Patient set ()
- Encounter set ()
- Previous query ()

7.2 Adding a folder and items

7.2.1 Folders and Items Draggable to Workplace


Folders and Items can be dragged and dropped (copied) into the 'Workplace' panel from other panels in the i2b2 Web Client.

A Concept Folder() and all its concepts (children), a Modifier Folder (), a concept (), a modifier () can be dragged to the 'Workplace' from the 'Navigate/Find Terms' panel.

A Patient set (), An Encounter set (), A previous query () can be added to the 'Workplace' by dragging from the 'Previous Queries' panel.

7.2.2 Adding a Folder


1. Using the *right mouse button*, click on your user folder in the 'Workplace'.
2. Select **New Folder and enter a name. Click Ok**
3. The new folder will appear in the 'Workplace' panel.

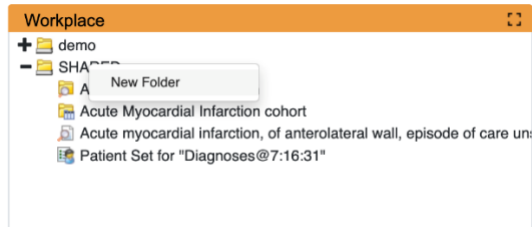

7.2.3 Adding Items

7.2.3.1 Adding Concept, Concept Folder, Modifier, and Modifier Folder

Adding from 'Terms' tab


1. In the Find Terms in 'Navigate/Find Terms' component, search for the item(s) to add.
2. Highlight the concept folder() or other type (concept, modifier, modifier folder) of item by clicking on its name.
3. While holding the *left mouse button* down, drag the item over to the 'Workplace' component.
4. Drop the item into either your workplace folder or the shared folder.
5. The item will now display in the 'Workplace' panel with the same name.

7.2.3.2 Adding Previous Query, Patient Set, and Encounter Set



'Queries' Panel

The **'Queries' panel** displays the queries that were run by the user. Each previous query contains a set of results; within these results are the result options that were selected when the query was originally run. The result options include Patient set, Encounter set, Number of patients as well as various other patient breakdowns

Adding Previous Query from 'Queries' Panel

1. Highlight the **previous query** () by clicking on the name of the query.
2. While holding the *left mouse button* down, drag the item to the **'Workplace' panel**.
3. Drop the previous query into either your workplace folder or the shared folder.
4. The previous query will now appear in the workplace with the same name.

Adding Patient Set or Encounter Set from 'Queries' Panel

1. Click on the plus sign (**+**) next to the name of the previous query that contains the patient set.
2. Click on the plus sign (**+**) next to the results folder.
3. Highlight the **Patient set** () or the **Encounter set** () by clicking on its name.
4. While holding the *left mouse button* down, drag the **Patient Set** or the **Encounter set** to the '**Workplace**' panel .
5. Drop the patient or encounter set into either your workplace folder or the shared folder.
6. The patient set will now appear in the workplace.

7.3 Navigation

7.3.1. Expand / Collapse a Folder


Users can expand or collapse a folder by clicking on the plus sign (**+**) or negative sign (**-**) next to the folder icon.

Expand (Open) a Folder

1. Click on the **plus sign** (**+**) next to the folder in the '*Workplace*' panel.
2. The folder will open and all the items in the folder will display.

Collapse (Close) a Folder

1. Click on the **negative sign** (**-**) next to the folder in the '*Workplace*' panel.
2. The folder will close and all the items in the folder will display.

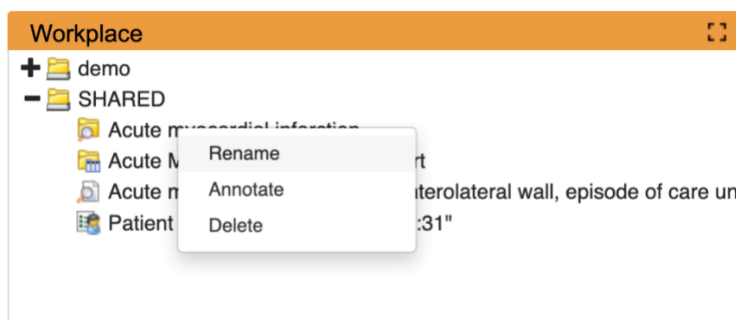

7.3.2 Rename a Folder or an Item

Users can easily **rename** an existing folder or an item.

1. Using the *right mouse button*, click on the folder or item you want to rename. Select **Rename**.
2. In the text box change the current name to the new name. Click **Ok**

3. The folder will appear in the 'Workplace' panel with the new name.

7.3.3 Annotate a Folder or item

Users can change the default **annotation** for any of their folders or items in the 'Workplace' panel.

1. Using the *right mouse button*, click on the folder or item you want to annotate. Select *Annotate*.
2. In the text box enter the text you want to appear when the tooltip is displayed.
6. Click on the **OK button**.
7. Using the *mouse*, hover over the folder or item you just annotated.
8. The text you entered will display as a tool tip.

7.3.4 Delete a Folder or item

Users can remove a folder from the 'Workplace' panel.

WARNING:


All items in the folder will be deleted when the folder is removed. It is recommended that you move those items you want to keep to another folder.

1. Using the *right mouse button*, click on the folder or item you want to delete. Select *Delete*.
2. Click on the **OK button to confirm**.
3. The folder and all of its contents will no longer appear in the 'Workplace' panel.

7.4 Using Items in 'Workplace' Panel






7.4.1 Items Draggable from Workplace


Items can be dragged and dropped (copied) from the 'Workplace' panel to other panels in the i2b2 Web client.

- Ontology Folder () and all its concepts (children)

NOTE:

User folders cannot be dragged from the 'Workplace' panel.

- Modifier Folder () and all its modifiers (children)
- Concept ()
- Modifier ()
- Previous Query ()
- Patient Set ()

- Encounter Set ()