

1.6.08 Release Notes

Release Notes for i2b2 Version 1.6.xx

These release notes address resolutions for issues found and information about the i2b2 version 1.6.00 to 1.6.08. Unless otherwise specified, these notes apply to the following operating systems:

- Apple Macintosh OS x 10.4/5
- Apple Macintosh OS x 10.6 (Snow Leopard)
- Microsoft Windows XP Professional
- Microsoft Windows 7
- Microsoft Windows Server 2003

Additional Information

Included in this document are the following sections:

Section	Includes
Installation	Information related to installing the release.
Third Party Software	Software that is utilized by i2b2 but is owned, maintained and licensed by other companies is referred to as "Third Party Software". Any version changes made or needed to be made to the users system will be noted in this section.
New Features	A list of new enhancements will be included in this section. A brief description and a reference to the specification will be provided.
Updates	Any changes that are not considered an enhancement will be included in this section. This may be a fix for a problem found or it may be a change to existing code but not related to a problem.
Licensing	i2b2 licensing information
Trademarks	Trademark information

Installation

This release note applies to you if you are upgrading to 1.6.xx from a computer running an earlier version of i2b2. If you are installing 1.6.xx from a computer running the i2b2 Workbench, ensure that the workbench and all eclipse applications are closed prior to installing this release.

WHERE DO I START?

Prior to setting up the i2b2 server, workbench or VM Image it is highly recommended that you read one of the following documents located in the documentation folder.

1. Hive_Installation_Guide.pdf	use for new installations
1. Data_Installation_Guide.pdf	use for new installations of the demo and metadata

1. i2b2Workbench_Install_Guide.pdf	use for installing the i2b2 Workbench (Client)
1. vmware_Install_Guide.pdf	use for installing the VM Image

All of these documents will explain what needs to be done during the setup and in what order.

Release 1.6.08

Release 1.6.08 contains changes to the i2b2 Workbench (eclipse client), Web client, and server therefore it is necessary to update only those three components. There are no changes to the database scripts.

Third Party Software

No changes to third party software are included in release 1.6.08

New Features

There are no new features in Release 1.6.08.

Updates

Data:

There are no changes to the data or the database scripts.

Fixes:

The following fixes were added to this release.

Web Client: Modifiers not always displaying correctly in Query Tool

When a modifier was included in a query, the web client did not always display the name correctly in the panel of the Query Tool. Instead of displaying the name of the associated concept with the modifier in brackets it was displaying the name of the modifier folder / container for that modifier. The correct results were returned as the correct <item_key> was sent in the request xml.

STATUS/RESOLUTION: *This issue has been resolved. The name of the item now appears correctly when a modifier is dropped in a panel within the Query Tool.*

Web Client: Views did not always appear correctly when resizing panels

In the previous release we introduced a new enhancement that allowed users to resize the two i2b2 panels in the Web Client (see notes for Release 1.6.07). Occasionally when a user resized the new panels and then refreshed their browser window the views within the two panels would not be sized correctly. The issue was random and difficult to reproduce.

STATUS/RESOLUTION: *This issue has been resolved. The views will now appear correctly when a browser is refreshed after resizing the panels. Every time a user logs into the Web Client the panels and views will display in the default setting; the left panel will be 40% of the page and the right panel will be 60%.*

CRC: Error returned when name contains a square bracket

The CRC would return an error when running a query that contained a concept or modifier with a square bracket in its path name. The right and left square brackets are special characters and the CRC would not handle them correctly when the SQL statement for the query was generated.

STATUS/RESOLUTION: *This issue has been resolved. All special characters will now be "escaped" in the SQL statement generated when running a query.*

ONT: Error returned when name contains a square bracket

Similar to the same issue with the CRC the Ontology would return an error when a square bracket was in the name or path of a concept or modifier. This was most notable when during a search in Find Terms.

STATUS/RESOLUTION: *This issue has been resolved. All special characters will now be "escaped" in the SQL statement when a request is received by the Ontology.*

Workbench: Date format sent in request message needs time zone

The date format sent in the XML request messages does not have the time zone information. This affects queries that are run with a date constraint and when generating a timeline with a date constraint (PDO request).

STATUS/RESOLUTION: *This issue has been resolved. The date format sent for the date constraints has been updated to include the timezone*

EXAMPLE: Correct date format in xml for date constraint

```
<constrain_by_date>
<date_from>1999-12-30T00:00:00-05:00</date_from>
<date_to>2005-01-01T00:00:00-05:00</date_to>
</constrain_by_date>
```

CRC: Temp table used for "simple queries"

On occasion a temp table was being used when running a "simple query". In this instance, a simple query is one in which the only result type requested and returned is the number of patients (patient count). The use of a temp table in this scenario is less efficient than a simple select statement. It is important to note that this issue only caused a problem with the efficiency of running a query and did not affect the results returned.

STATUS/RESOLUTION: *This issue has been resolved. A temp table will no longer be used when running a simple query for the number of patients that meet the criteria defined.*

Release 1.6.07

Release 1.6.07 contains changes to the i2b2 Web client; therefore it is necessary to update only the Web client. There are no changes to the i2b2 Workbench, server or the database.


Third Party Software

No changes to third party software are included in release 1.6.07

New Features

WEB CLIENT: Ability to "Resize" panels

In release 1.6.07, users can now resize the two panels in the i2b2 Web client. This is accomplished by dragging the vertical bar that now appears in between the two panels to either the right or the left. When you drag it to the right, the Navigate Terms, Find Terms, Workplace, and Previous Query views will become wider. When you drag it to the left, the Query Tool and Query Status views (or plugin viewer) will become wider.

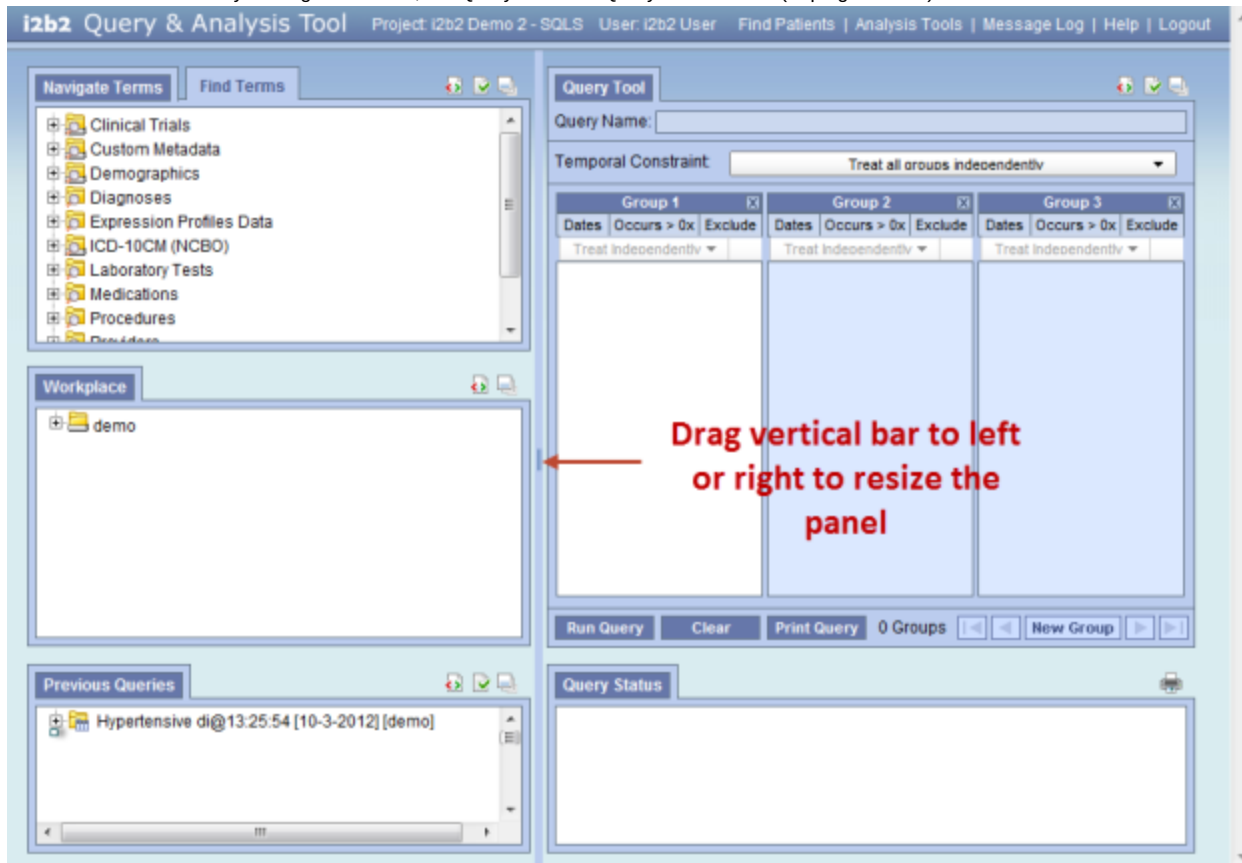

WEB CLIENT: Edits made to a query will clear the results in the Query Status view

To prevent confusion and possible misinterpretation the results in the Query Status view will be cleared whenever a change is made to a query in the Query Tool. For example, when you drag a previous query to the Query Tool the results for that query will show in the Query Status. If you add a concept or change a constraint then the results in the query status no longer accurately represent the query defined in the Query Tool. Therefore the results are cleared from the Query Status view and the name of the previous query will no longer appear at the Query Name in the Query Tool.

Updates

Data:

There are no changes to the data or the database.

Fixes:

The following fixes were added to this release.

Web Client: Queries with a modifier did not always send the correct key in the XML message

When a modifier was included in a query, the web client did not always send the correct <item_key> in the XML request message. The problem would only occur if the modifier was a child of another modifier (folder). In this instance the <item_key> for the modifier folder was being sent instead of the key for its associated concept, which resulted in zero patients being returned.

STATUS/RESOLUTION: *This issue has been resolved. The correct <item_key> is now sent regardless of whether or not the modifier is a child of another modifier.*

Release 1.6.06

Release 1.6.06 contains changes to the i2b2 Web client and server; therefore it is necessary to update only those two components. There are no changes to the i2b2 Workbench or the database.

Third Party Software

No changes to third party software are included in release 1.6.06

New Features

WEB CLIENT: "By Value" now defaults in value constraint window

A simple change was made to the value constraint window in that prior to 1.6.06 the selection of "No Value" would default when the value constraint window opened. Now in 1.6.06 the default is "By Value" or "Search within text".

WEB CLIENT: Bar graph added to value constraint window

In release 1.6.06 the value constraint window for numerical terms has been modified to display a bar graph with the range of available values (see screenshot shown below).

- *"numerical term" is defined as a term that has the datatype of Integer, PosInteger, Float, PosFloat defined as the DataType in the c_metadataxml column in the i2b2 table.*

Now there are two ways in which a user can enter a value:

1. Manually enter it at "Please enter a value"
2. Click on the bar graph to enter the value for that segment.

If you click on the bar graph the operator will automatically update to either "Less Than or Equal To" or "Greater Than or Equal to". The operator that appears depends on the segment you clicked on.

Example: The green segment is considered the "Normal Range".

- If you click on a segment below the normal range then the operator will change to "Less Than or Equal To".
- If you click on a segment above the normal range then the operator will change to "Greater Than or Equal To".
- *All of the operators are still available and can still be selected by clicking on the down-arrow to display a drop-down list. The process of updating the operator automatically is designed to save the user key strokes by having the most common selection default.*

Choose value of CPK (Test:CPK)

Searches by Lab values can be constrained by the high/low flag set by the performing laboratory, or by the values themselves.

☐ No Value

☐ By FLAG

☒ By Value

Please select operator:
LESS THAN (<)

Please enter a value:

Click on a bar segment to help specify a value or range:
Range in u/l

0 61 150 400

Units = u/l

OK Cancel

In order for the bar graph to appear the concept or modifier has to have metadatatxml entered in the metadata table for that item. To be more specific, the ranges have to be defined in the following xml tags:

```
<LowofLowValue></LowofLowValue>
<HighofLowValue></HighofLowValue>
<LowofHighValue></LowofHighValue>
<HighofHighValue></HighofHighValue>
<LowofToxicValue></LowofToxicValue>
<HighofToxicValue></HighofToxicValue>
```

The above tags are not new but in order for the bar graph to appear it is necessary that values be entered. The metadatatxml for the above example (screenshot) is as follows:

```
<?xml version="1.0"?>
<ValueMetadata>
<Version>3.02</Version>
```

```

<CreationDateTime>04/15/2007 01:22:54</CreationDateTime>
<TestID>CPK</TestID>
<TestName>CPK</TestName>
<DataType>PosInteger</DataType>
<CodeType>GRP</CodeType>
<Loinc>2157-6</Loinc>
<Flagstouse>HL</Flagstouse>
<Oktousevalues>Y</Oktousevalues>
<MaxStringLength></MaxStringLength>
<LowofLowValue>0</LowofLowValue>
<HighofLowValue>61</HighofLowValue>
<LowofHighValue>150</LowofHighValue>
<HighofHighValue>400</HighofHighValue>
<LowofToxicValue></LowofToxicValue>
<HighofToxicValue></HighofToxicValue>
<EnumValues></EnumValues>
<CommentsDeterminingExclusion>
<Com></Com>
</CommentsDeterminingExclusion>
<UnitValues>
<NormalUnits>u/l</NormalUnits>
<EqualUnits>u/l</EqualUnits>
<ExcludingUnits></ExcludingUnits>
<ConvertingUnits>
<Units></Units>
<MultiplyingFactor></MultiplyingFactor>
</ConvertingUnits>
</UnitValues>
<Analysis>
<Enums />
<Counts />
<New />
</Analysis>
</ValueMetadata>

```

You will notice that the ranges that appear in the screen shot correspond to the values entered at the "...Value" tags.

Updates

Data:

There are no changes to the data or the database.

Fixes:

The following fixes were added to this release.

ONT: Hidden and Synonym modifiers are appearing in Navigate Terms view.

Both hidden and synonym modifiers are appearing in the Navigate Terms view even though neither option is selected.

STATUS/RESOLUTION: *This issue has been resolved. Both hidden and synonym terms will only appear when they are selected in the Options dialog.*

WEB CLIENT: Synonyms are not always appearing in Find Terms view.

In the sent XML message sometimes the web client was sending synonyms="false" even though the option to "Show Synonymous Terms" was selected. As a result synonyms were not displaying in the Find Terms view.

STATUS/RESOLUTION: *This issue has been resolved. The correct setting will be sent in the xml message to the server. Synonyms will now display in the Find Terms view as long as the option to show them is selected.*

WEB CLIENT: Panel Timing is not always displayed correctly.

In the Web client if there are more than three panels being used in the Query Tool AND different panel timings are defined then the correct panel timing may not display when you scroll through the panels. The timing appears to shift to the next panel. This is only a display issue because when the query is run the correct panel timing is sent in the xml message.

STATUS/RESOLUTION: *This issue has been resolved. The selected panel timing will always display with the appropriate panel regardless of how many panels are in use.*

Release 1.6.05

Release 1.6.05 contains changes to the i2b2 Workbench (eclipse client), Web client, server, and database scripts therefore it is necessary to update all four components.

Third Party Software

No changes to third party software are included in release 1.6.05

New Features

Workbench: "Number of patients" is checked by default

In release 1.6.05 the "Number of patients" analysis type is now checked by default in the query tool.

Workbench: "Timeline" checked by default if user has the appropriate permissions to view timeline data

Prior to release 1.6.05, the Analysis Type of "Timeline" was checked off by default when a user ran a query. Since the availability of the Timeline view is restricted to a user's access level (data_ids or higher) this default was not appropriate for all users. Therefore in release 1.6.05 the Analysis Types that default when running a query has been modified to look at the user's role.

Users Role	Default Analysis Type
DATA_LDS, DATA_DEID, DATA_PROT	Number of Patients and Timeline
DATA_AGG, DATA_OBFSC	Number of Patients

- *This change only affects the workbench. The web client will continue to default to Number of Patients.*

SQL Server Indexes Optimized

The SQL Server indexes were updated for optimization during searches and querying.

Updates

Data

The only changes made to the database scripts were to the stored procedures. The changes were made as part of the fixes for the issue with the CRC uploader ("CRC: Error returned when importing items into Oracle database").

Fixes:

The following fixes were added to this release.

WEB CLIENT: Not able to expand project users in the Admin Tool.

In the Admin Tool, the list of users would not display if you used the plus sign to expand the list. In addition, the admin tool would freeze and users had to log back into it to go to any other pages. The problem only occurs when using the plus / minus sign next to the name it will not occur if you click on the name.

STATUS/RESOLUTION: *This issue has been resolved. List of project users will expand if click on the plus sign or the name.*

WEB CLIENT: Queries with occurrence greater than 1 may return an error.

In the Web Client, an error is returned if running a query with an occurrence constraint defined. This is only happening if the Project Request plug-in is removed from the Web Client. An error will not be returned if the plug-in remains installed.

STATUS/RESOLUTION: *This issue has been resolved. Queries with an occurrence constraint will work correctly regardless of whether or not the Project Request plug-in is installed.*

CRC: Error returned when importing items into Oracle database.

The CRC is returning an error when trying to import files into an Oracle database using the import view.

STATUS/RESOLUTION: *This issue has been resolved. Observation_facts, patient_mapping and event_mapping files can now be imported into an oracle database.*

ONT: SQL server only – Indexes Optimized.

The correct indexes were not always used by the ontology cell without specification in JDBC Source statement. This is only a problem in SQL Server and is not a problem in Oracle.

STATUS/RESOLUTION: *This issue has been resolved. The ont-ds.xml file was updated and the correct indexes are now being used.*

Release 1.6.04

Release 1.6.04 contains changes to the i2b2 Workbench (eclipse client), Web client and server therefore it is necessary to update all three components. Minor edits were made to the data in the i2b2 metadata. No modifications were made to the database structure.

Third Party Software

No changes to third party software are included in release 1.6.04

New Features

i2b2 release 1.6.04 does not contain any new features.

Updates

Data

A few of the items in the i2b2 metadata table were missing values in the M_EXCLUSION_CD column and this resulted in the concepts not being copied to the concept_dimension table when a synchronization in the edit terms view was performed.

No changes were made to the database structure.

Fixes:

The following fixes were added to this release.

Queries with occurrence and exclusion constraints not always working correctly

Queries that contain at least one panel that has an occurrence constraint and is excluded will not return the correct number of patients. This is only an issue if both constraints appear on the same panel.

STATUS/RESOLUTION: *This issue has been resolved. Excluded panels that contain an occurrence constraint will now return the correct number of patients.*

Query-in-Query produces a "union" rather than a "join" when running the query.

A query-in-query is a query that is run using a previous query in at least one panel. There is an issue if more than one previous query is used in at least two panels (groups). It appears that when the query-in-query is run, a union of the two panels occurs rather than a join. As a result the correct number of patients is not returned.

STATUS/RESOLUTION: *This issue has been resolved. Multiple previous queries in multiple panels will now be joined. The correct number of patients will be returned.*

Query: Square bracket ([']) in an item's name will return an error

Running a query that contains at least one item that has a square bracket ([']) in its name will result in an error from the server.

STATUS/RESOLUTION: *This issue has been resolved. Square brackets ([']) in the name will no longer return an error.*

Edit Tool: Synchronize process does not run

The system appears to "hang" when running the Synchronize All or Update process in the Edit Terms view.

STATUS/RESOLUTION: *This issue has been resolved. The synchronization process runs correctly.*

Token for User Session prematurely times out

When a user logs into the i2b2 Workbench or Web client a token for that session is assigned. There is a problem in which this token is timing out too early and as a result the user has to log out and log back into the workbench / web client.

STATUS/RESOLUTION: *This issue has been resolved. The session token will now time out in the designated time frame.*

Release 1.6.03

Release 1.6.03 contains changes to the i2b2 Workbench (eclipse client), Web client and server therefore it is necessary to update all three components. It is also important to note that additional data has been added to the demo and metadata to support the new features.

Third Party Software

No changes to third party software are included in release 1.6.03

New Features

i2b2 release 1.6.03 contains new features that have been outlined below.

Large Text Searches

In the previous release, 1.6.02 we introduced a new feature called "String Constraints". Basically a string search ("constraint") is similar to other value constraints in that users can narrow down their search criteria based on a small string of text stored in the tval_char column of the observation_fact table. (For more information on string constraints see the notes for Release 1.6.02).

In release 1.6.03 we now introduce the ability to narrow down your search based on text that is stored in the observation_blob column of the observation_fact table. This new feature is called "Large Text Searches" as the data that is stored in this column is typically a large blob of text such as a report.

To the user the process defining a text search, whether it is a "small string" or "large text string" is pretty straight forward and similar. There are a few minor differences they will see. For example:

- Small string search users can define the value operator (example; contains, starts with, ends with, and exact).
- Large text string search works similar to other search engines; Specifically users can define database search operators in the string of text they enter.

The metadata.xml defined for the concept or modifier will determine which value constraint window will display to the user. The metadata.xml will also determine the <value_type> that is sent in the xml request message. This <value_type> will define whether or not the tval_char or the observation_blob column is searched for the defined text.

IMPORTANT: If you are upgrading from 1.6.02 to 1.6.03 you need to run the following against your database.

```
insert into QT_PRIVILEGE(PROTECTION_LABEL_CD,DATAPROT_CD,HIVEMGMT_CD) values  
( 'SETFINDER_QRY_WITH_LGTEXT','DATA_DEID','USER');
```

WEB CLIENT: Auto-generate "Timeline"

In the past, when running a query with the result type of Timeline selected it would simply return a patient set. In 1.6.03 this feature has been updated so that now when it is selected, the timeline plug-in will open and automatically display the results in the timeline. This feature is already available in the i2b2 Workbench (eclipse client) so this change is to add it to the web client.

- ***Currently modifiers are not supported in the timeline plug-in for the web client.***

WEB CLIENT: View text (reports) in the timeline

In 1.6.03 users can now view text (reports) from within the timeline. When double clicking on the tick mark that appears for the report a new window will open displaying the text. This feature is already available in the i2b2 Workbench (eclipse client) so this change is to add it to the web client.

- ***Currently encrypted notes are not supported and cannot be viewed in the timeline plug-in for the web client. They can be viewed in the i2b2 workbench (eclipse client) with the appropriate encryption key.***

WEB CLIENT: Print Query added to Query Status

The option to print a query has been added to the Query Status view of the Web Client. A new icon () will now appear in the top right corner of the Query Status view. A new window will open when you click on this icon. In this window all the information regarding the query and its constraints will appear in the header of the document. The results returned will appear below the header in the body of the report.

Updates

Data

In release 1.6.03, no changes were made to the structure of the database. The only change was to add examples to the metadata, concept_dimension, modifier_dimension, and observation_fact tables to support the new large text constraint.

In addition to adding new data examples the metadataxml was updated with a new data type of "largestring". An example of this new metadataxml is as follows:

```
<?xml version="1.0" encoding="UTF-8"?>
<ValueMetadata>
<Version>3.2</Version>
<CreationDateTime>2011-12-19T13:32:16.198-04:00</CreationDateTime>
<TestID>LCS-I2B2:XR_RPT_DID</TestID>
<TestName>De-Identified X-Ray Reports</TestName>
<DataType>largestring</DataType>
<MaxStringLength/>
<Flagstouse/>
<Oktovalues/>
</ValueMetadata>
```

Fixes:

The following fixes were added to this release.

Queries with exclusion not always working correctly

Queries that contain multiple items in a panel defined to be excluded will not return the correct number of patients. This appears to be happening due to the order in which the inversion is taking place when the SQL for the query is generated and run.

STATUS/RESOLUTION: *This issue has been resolved. Excluded panels will always be performed at the end of the SQL statement. All items in the panel will now be excluded.*

Hidden terms not returned in searches

Hidden Terms were not being returned in searches even though it was selected in the options.

STATUS/RESOLUTION: *This issue has been resolved. Hidden terms will now be returned when "Show hidden terms" is selected in the options dialog.*

Error returned if an item in a previous query contained a Hidden Term

Using a previous query that contains a term that has been changed to be a hidden term will return an error. This is only an issue if the term was active when the query was originally run and later changed to hidden.

STATUS/RESOLUTION: *This issue has been resolved. Previous queries with a hidden term can be dropped in the query tool and the new query will run..*

Find Terms: First term in list sometimes inherits color of terms in the previous search

Terms will appear in blue if they are a synonym and red if they are a hidden term. If you search for a term that has a term highlighted in red or blue then when you do the next search the first term in the list will inherit the color from the previous search.

STATUS/RESOLUTION: *This issue has been resolved. All terms now appear in the correct color regardless of the results from previous searches.*

Macintosh: Expand / Collapse did not always work in Previous Query

In the Previous Query view, users cannot expand or collapse folders once you highlighted an item in the view. You can still open and close a folder by double clicking on the name but you cannot do it by clicking on the icon (arrow) next to the name.

STATUS/RESOLUTION: *This issue has been resolved. Users can open and close a folder by clicking on the icon or double clicking on the name*

- *This was only an issue when using a Macintosh to log into the i2b2 Workbench. The process worked correctly if using Windows.*

Release 1.6.02

Release 1.6.02 contains web client, server and oracle database changes so it is important to update all that apply to you. This release does not contain any changes to the i2b2 Workbench.

Third Party Software

No changes to third party software are included in release 1.6.02

New Features

i2b2 release 1.6.02 contains new features that have been outlined below.

String Constraint

In the 1.6.02 release there is now an option to constrain a query by a string of text. Similar to other text and numeric constraints, a user can narrow down their search criteria when setting up the query in the query tool view. The main difference between this new constraint and the existing ones is the string constraint allows users to enter a string of text as opposed to an enumerated or numerical value.

The operators for a string constraint are defined below.

Operator	Description
Starting with	The results returned will contain only those facts whose value begins with the text entered in the search box.
Ending with	The results returned will contain only those facts whose value ends with the text entered in the search box.
Containing	The results returned will be those facts whose value contains the text entered in the search box. • <i>This constraint will return the broadest range of results.</i>

Exact	The results returned will contain only those facts whose value matches exactly to the text entered in the search box. • <i>This constraint will return the narrowest range of results.</i>
--------------	--

Defining the string constraint works the same way as other text and numeric constraints. Once you drag the concept or modifier to a panel in the query tool view the appropriate constraint window will open. The only requirement is that the concept or modifier has to have metadataxml entered in the metadata table for that item. This requirement has not changed. The only thing that has changed is a new data type has been added to the metadataxml. Information about this new data type can be found in the Data section of these release notes.

CDATA marker added to the xml message

The 1.6.02 release, the CDATA marker has been added to the xml messages being sent from the i2b2 Web Client and the Workbench. This new section will appear for those concepts and modifiers that can have a value constraint defined when running a query or generating a timeline. An example of the change is shown below.

```
<item>
<hlevel>3</hlevel>
<item_name>Social history [Contains: smokes]</item_name>
<item_key>\\i2b2_REP\\i2b2\\Reports\\Medical History\\Social History{color:#0000ff}</item_key>
<tooltip>i2b2 \ Reports \ Medical History \ Social History</tooltip>
<class>ENC</class>
<item_icon>LA</item_icon>
<item_is_synonym>>false</item_is_synonym>
<constrain_by_value>
<value_type>TEXT</value_type>
<value_operator>LIKE[contains]</value_operator>
<value_constraint><![CDATA[smokes]]></value_constraint>
</constrain_by_value>
</item>
```

Lockout / Obfuscation Parameters

In the 1.6.02 new parameters were added to the crc.properties file to support the lockout process for obfuscated users as well as the obfuscation process itself.

In addition to the new parameters a modification was made to one of the existing parameters (*edu.harvard.i2b2.crc.lockout.setfinderquery.count*). The value of -1 is now supported, which means if this parameter is set to -1 then obfuscated users will not be locked out.

This modification as well as the new parameters is outlined below. More detailed information can be found on the i2b2 Community Wiki.

Parameter	Description
edu.harvard.i2b2.crc.lockout.setfinderquery.count	The number of times a user is allowed to run a query with the same results before they are locked out. <i>CHANGE: The value of -1 is now supported.</i>
edu.harvard.i2b2.crc.lockout.setfinderquery.zero.count	Defines whether or not queries with a set size of zero will be counted when determining whether or not to lockout a user.
edu.harvard.i2b2.crc.setfinderquery.obfuscation.count.sigma	The standard deviation to be applied when displaying the obfuscated results (set size).
edu.harvard.i2b2.crc.setfinderquery.obfuscation.breakdowncount.sigma	The standard deviation to be applied to the set size for the "patient breakdowns" analysis types.

edu.harvard.i2b2.crc.setfinderquery.obfuscation.minimum.value	The threshold (minimum value) a result must meet. Any value below this number will return zero as the new obfuscated set size.
---	--

Updates

Data

In release 1.6.02, no changes were made to the structure of the database. The only change was to add examples to the metadata, concept_dimension, modifier_dimension, and observation_fact tables to support the new string constraint.

In addition to adding new data examples the metadata.xml was updated with a new data type of "String". An example of this new metadata.xml is as follows:

```
<?xml version="1.0" encoding="UTF-8"?>
<ValueMetadata>
<Version>3.2</Version>
<CreationDateTime>2011-10-19T13:32:16.198-04:00</CreationDateTime>
<TestID>MED:SIG</TestID>
<TestName>Rx Sig</TestName>
<DataType>String</DataType>
<MaxStringLength>50</MaxStringLength>
<Flagstouse/>
<Oktousevalues/>
</ValueMetadata>
```

Fixes:

The following fixes were added to this release.

Memory error: not able to load PDO

When running a query to obtain a patient list an out of memory error will occur if the database contains a large amount of data. In addition, the patient data object (PDO) is never fully created and will not open in the previous query view.

STATUS/RESOLUTION: *This problem was introduced with the enhancement to allow for custom columns in the PDO and only affects large databases. This issue has been resolved. Creating and loading the PDO will no longer cause an out of memory error.*

Not able to search for terms with an '&' in the name

In the i2b2 Web client users are not able to search for terms with an ampersand '&' in the name. The term would never be returned with the results.

STATUS/RESOLUTION: *This issue has been resolved. Terms with an '&' in the name are now returned in the Find Terms View.*

Demo Database Scripts: Passwords inconsistent

In the scripts for the i2b2 demo database the passwords for the user setup during the install process are inconsistent.

STATUS/RESOLUTION: *This issue has been resolved. All of the above users now have the same password ('demouser').*

- *The password is only changed for the i2b2 demo database. Your individual administration and user passwords for your own databases will not change.*

PM: Project-User parameters are not saved with the project

The project-user parameters entered into Admin are not being saved with the project id. Therefore they are being associated with all projects instead of the one intended.

STATUS/RESOLUTION: *This issue has been resolved. Project-User parameters are now saved with the appropriate project id and are specific to the project and the user.*

- *If you have entered project-user parameters in version 1.6 or 1.6.01 then you will need to fix them after you upgrade to 1.6.02. Once you upgrade to 1.6.02 the parameters with 'undefined' as the project id will no longer appear in Admin. Therefore you will need to go into the PM_PROJECT_USER_PARAMS table and update the PROJECT_ID column with the correct project id.*

Release 1.6.01

Release 1.6.01 contains web client, server and oracle database changes so it is important to update all that apply to you. This release does not contain any changes to the i2b2 Workbench.

Third Party Software

No changes to third party software are included in release 1.6.01

New Features

i2b2 release 1.6.01 does not contain any new features.

Updates

Data

In release 1.6.01, no changes were made to the structure of the database. The only change was made to the metadata tables for the oracle database, which have been updated to remove the values in the C_TOTALNUM column. The SQL version was correct and therefore no changes are made to the SQL version.

Fixes:

The following fixes were added to this release.

Ontology: Need to escape apostrophe when searching for a term

Terms that contain an apostrophe are not returned when performing a search in the Find Terms view.

STATUS/RESOLUTION: *This issue has been resolved. The apostrophe in the term name or code will now be escaped and all the appropriate terms will be returned during a search.*

CRC: Need to escape apostrophe in SQL statements

Zero patients are returned when performing a query using at least one term that contains an apostrophe in the path name.

STATUS/RESOLUTION: *This issue has been resolved. All apostrophes will be escaped in SQL statements and the correct number of patients are returned when running a query.*

Project Management: Not able to add new users

An error is returned when trying to add a new user in the Project Management application.

STATUS/RESOLUTION: *This issue has been resolved. New users can be added to the PM.*

Web Client: The "Help" contains references to yui instead of local

The links in the on-line help contain references to the external yui. These links should be defined as local.

STATUS/RESOLUTION: *This issue has been resolved. The references in the on-line help file have been changed to local.*

Web Client: Scrollbar is missing from some analysis plug-ins

The horizontal and vertical scrollbars are missing from the following analysis plug-ins: Project Requests, Demographics (1 Patient Set), and Demographics (2 Patient Sets).

STATUS/RESOLUTION: *This issue has been resolved. The scrollbars have been added to the analysis plug-ins.*

Web Client: Analysis tools expand into Ontology views

When select Analysis tools from the menu bar, the plug-in viewer overlaps into the ontology views. As a result the task icons (options, xml, resize workspace) move over so that they appear above the plug-in viewer instead of the ontology views.

STATUS/RESOLUTION: *This issue has been resolved. The plug-in viewer no longer overlaps into the ontology views.*

Web Client: Not able to open a modifier container or folder

Modifier containers and folders would not expand to display the modifiers contained within them.

STATUS/RESOLUTION: *This issue has been resolved. All containers and folders will now open.*

Release 1.6.00

Release 1.6.00 contains changes to the i2b2 Workbench (eclipse client), Web client and server so it is necessary to update all three components.

Third Party Software

The YUI, which is used by the i2b2 Web Client, was upgraded to version 2.8.2

Information about YUI can be found on the YUI Library website (<http://developer.yahoo.com/yui>)

New Features

i2b2 release 1.6.0030 contains new features that have been outlined below.

Temporal (Query) and Panel Timing Constraints

In the 1.6 there is an option to **query by visit**. In the visit table of the data repository, there is a listing of visits that are tied to patients and facts in the observation_fact table through the encounter_num column. Both the temporal (query) and panel constraints are very similar in that both utilize the visit / encounter information along with the concepts defined to determine if a patient meets the search criteria.

The main difference between a **temporal** and a **panel timing constraint** is that temporal constraints are defined at the *query level* and panel timing constraints are defined at the *panel level*.

- When a query is requested to run and it is constrained by the **"Same Financial Encounter"**, the data will be tied together by the *encounter_num* and positive "hits" will only yield patients who have these observation_facts occurring together in one of their visits.
- When a query is requested to run and it is constrained by the **"Instance will be the same"**, the data will be tied together by the *encounter_num* and instance_num; positive "hits" will only yield patients who have these observation_facts occurring together in same instance of their visit.
- *In addition to the encounter_num and instance_num matching the following additional columns have to match in order to be considered the "same instance" patient_num, concept_cd, start_date, and provider_id*

This similar logic is applied to the panel timing. The main difference being that the timing is defined for each individual panel as opposed to the entire query. To avoid confusion the following outlines the labels seen in the client and the associated timing of any or same.

Timing Level (Query or Panel)	Label in Client	Timing (ANY or SAME)
*** Query Timing	Treat all groups independently	ANY
*** Query Timing	Selected groups occur in the same financial encounter	SAME
Panel Timing	Treat Independently	ANY
Panel Timing	Occurs in Same Encounter	SAME

*** In the web client "Query Timing" is called "Temporal Constraint". Although they are labeled differently the functionality is the same.

Example 1:

Patient has a medication with a dose greater than 500 mg **and** a medication that is taken once a day documented at **any time** in their history.

Query Tool

Query Name:

Temporal Constraint: Treat all groups independently

Group 1			Group 2			Group 3		
Dates	Occurs > 0x	Exclude	Dates	Occurs > 0x	Exclude	Dates	Occurs > 0x	Exclude
Treat Independently			Treat Independently			Treat Independently		
s [Medication Dose > 500 mg]			Frequency = ("Once per day")					

2 Groups

RESULTS: there are 47 patients who have had a medication with a dose greater than 500 mg and a frequency of once a day entered as a fact in the observation_fact table, where the medication dose and frequency have the same patient_num.

Example 2:

Patient has a diagnosis of Acute MI documented at *any time* in their history and they have a medication with a dose greater than 500 mg *and* a medication that is taken once a day documented during the *same visit*.

Query Tool

Query Name:

Temporal Constraint: Selected groups occur in the same financial encounter

Group 1			Group 2			Group 3		
Dates	Occurs > 0x	Exclude	Dates	Occurs > 0x	Exclude	Dates	Occurs > 0x	Exclude
Treat Independently			Occurs in Same Encounter			Occurs in Same Encounter		
Acute myocardial infarction			s [Medication Dose > 500 mg]			Frequency = ("Once per day")		

3 Groups

RESULTS: there are 4 patients who have had a medication with a dose greater than 500 mg and a frequency of once a day entered as a fact in the observation_fact table, where the medication dose and frequency have the same encounter_num AND the diagnosis of acute MI, the medication dose and frequency all have the same patient_num.

Example 3:

Patient has a diagnosis of Acute MI documented at *any time* in their history and they have a medication with a dose greater than 500 mg *and* a medication that is taken once a day documented during the *same instance*.

The screenshot shows the 'Query Tool' window. At the top, there's a 'Query Name' field. Below it is a 'Temporal Constraint' dropdown set to 'Items Instance will be the same'. The main area contains three groups, each with a header bar (Group 1, Group 2, Group 3) and a sub-header bar with 'Dates', 'Occurs > 0x', and 'Exclude' options. Group 1 has a dropdown set to 'Treat Independently' and contains the criteria 'Acute myocardial infarction'. Group 2 has a dropdown set to 'Items Instance will be the same' and contains the criteria 's [Medication Dose > 500 mg]'. Group 3 has a dropdown set to 'Items Instance will be the same' and contains the criteria 'Frequency = ("Once per day")'. Each group's criteria is preceded by a green box labeled 'one or more of these'. The groups are connected by blue boxes labeled 'AND'. At the bottom, there are buttons for 'Run Query', 'Clear', 'Print Query', and a status indicator '3 Groups'. Navigation arrows and a 'New Group' button are also present.

RESULTS: there are 3 patients who have had a medication with a dose greater than 500 mg and a frequency of once a day entered as a fact in the observation_fact table, where the medication dose and frequency have the same encounter_num AND instance_num AND the diagnosis of acute MI, the medication dose and frequency all have the same patient_num.

Encounter Sets

In the 1.6 Release users are now able to run a query to return a set of encounters. These encounter sets are similar to patient sets only they contain a list of encounter numbers for the patients meeting the criteria defined when running a query.

Query-in-Query

The ability to use a previous query within a new query is now available in the 1.6. Release. To use this feature simply drag the previous query to one of the panels (groups) in the query tool. Add your additional search criteria and run the new query.

Query Using Patient and/or Encounter Sets

Similar to the query-in-query enhancement, users can now use a patient or encounter set from a previous query in a new query. This is done by dragging the patient or encounter set to a panel (group) in the query tool and then adding the additional search criteria. When you run the query it will use the patient /encounter set to query against.

Performance Tracking

The ability to track and display performance metrics when running a query has been added to 1.6 Release. The tracking information can be turned on by setting the following project parameter.

Parameter Name: PM_ENABLE_PROCESS_TIMING

Parameter Value: DEBUG

When the performance tracking is turned on (set to DEBUG) the information will appear in the Query Status view when you run the query.

- *This change is only available in the i2b2 Web Client.*

Update Term Usage Utility

In the 1.6 Release a new utility has been added to the i2b2 Workbench. This new utility enables users to easily run an update process which will update each term with the total number of times the term is used in the observation_fact table.

These totals can be seen in the Navigate Terms view by turning on the "Enable patient counts" in the Options dialog. This also enables performance tuning by the query engine.

- *This view is available in the workbench and the process can only be run by a user who has the role of "Editor" associated to them.*

Ontology View – Updates to Options Dialog

The existing Options dialog in the Navigate terms view has been modified slightly in 1.6 Release. The following items have been removed from the pop-up menu and added to the options dialog.

- Enable Patient Counts
- Use Short Tooltips
- Show Concept Codes in tooltips

Each of the above items works the same as they did when they were included in the pop-up menu.

Analysis Breakdowns Added to Web Client

In the 1.6 Release the following Analysis breakdowns have been added to the web client.

- Gender patient breakdown
- Vital Status patient breakdown
- Race patient breakdown
- Age patient breakdown

These breakdowns are available when running a query; the user can check the query result type(s) for each of the breakdowns. The results can be visible in the Query Status view when the query has completed running. They can also be seen when you drag the previous query to the Query Tool view, the breakdowns / analysis types for that query will appear in the query status view.

Modifiers

The ability to use modifiers when querying and viewing data from within the i2b2 workbench and web client has been implemented in 1.6 Release. The usage of modifiers allows a single fact to be customized with an unlimited number of codes. For instance,

- A medication can be modified with the dose, route, and frequency.
- A procedure can be modified with a CPT modifier(s) indicating a two surgeon team with anesthesia by anesthesiologist (62 and AA).
- A diagnosis involving malignant neoplasms can be modified with the TNM staging information.
- Allergies and other diagnoses can be modified with the severity.

- *The changes made to the database to support modifiers are listed in the Update section under Data*

Unit Conversion

In the 1.6 Release the ability to run queries using units other than the normal unit has been added. This feature is beneficial for those numerical facts that are not stored in normalized unit_cd.

The user can enable the unit conversion process by setting the following Project Parameter.

Parameter Name: CRC_ENABLE_UNITCD_CONVERSION

Parameter Value: ON

- *This unit conversion option will slow the query, so for better query performance load the numerical fact values in the normalized units and do not enable this option.*

When the above parameter is set to 'ON', the unit conversion of the nval_num is calculated using the concept's metadata xml defined in the ontology cell (<ConvertingUnits/> , <MultiplyingFactor/>). The unit conversion in the query is done before applying the value constraint.

Updates

Data

The demo data and metadata tables have been updated to support the new enhancements. This includes but is not limited to the following:

- Populating the encounter_mapping table
- Synchronizing all tables that store the encounter_num
- Modifier codes and instance numbers were updated in the observation_fact table
- New modifier_dimension table
- New columns added to the metadata tables
- New column added to both the patient_dimension and visit_dimension tables

Changes to support Modifiers

A number of changes have been made to the database to support the use of modifiers. A new table has been created, two new columns have been added to the metadata tables and three new values have been added to the C_VISUALATTRIBUTES column.

NEW TABLE

The **modifier_dimension** table is similar to the *concept_dimension* table except it stores the modifiers.

NEW COLUMNS

New Column	Description	Example
M_APPLIED_PATH	Contains the path of the concept that the modifier term applies.	In this example the modifier <i>Severity: Mild</i> will apply to the root concept of <i>'Diagnoses'</i> and all its descendents. Modifier Term = Severity: Mild M_APPLIED_PATH = \i2b2\Diagnoses%
M_EXCLUSION_CD	A flag to indicate whether or not a modifier term is not associated with a specific concept.	
'X'	Modifier is excluded	
NULL	No exclusions	
In this example the modifier <i>Severity: Mild</i> will not apply to the root concept of <i>'Neoplasms'</i> and all its descendents. Modifier Term = Severity: Mild M_APPLIED_PATH = \i2b2\Diagnoses\Neoplasms (140-239)% M_EXCLUSION_CD = X		

It is important to note that when excluding a modifier term the table will contain multiple entries for that modifier.

The following is a snippet from the i2b2 metadata table. In this example, the modifier term called *Severity: Mild* will be applied to the Diagnoses concept and all its children **EXCEPT** for the following two concepts and all their children. (1) Neoplasms and (2) Mental Disorders

- This snippet only includes those columns that are relevant to the example.*

C_FULLNAME	C_NAME	M_APPLIED_PATH	M_EXCLUSION_CD
\Severity\Mild\	Severity: Mild	\i2b2\Diagnoses%	NULL
\Severity\Mild\	Severity: Mild	\i2b2\Diagnoses\Neoplasms (140-239)%	X
\Severity\Mild\	Severity: Mild	\i2b2\Diagnoses\Mental Disorders (290-319)%	X

NEW VISUAL ATTRIBUTES

New Values	Description / Use	Example
OA	Modifier <i>Container</i>	TNM Staging
DA	Modifier <i>Folder</i>	TNM Staging, Stage Classification

RA	Modifier <i>Term</i>	Stage II: Cancer is localized, advanced within area
-----------	----------------------	---

Documentation:

The 1.6.0030 release documentation can be found in a zip file located on the i2b2 Software page. This file includes documentation for the following:

- The architecture, design, installation and messaging documents for the CRC, Ontology, File Repository, PM and Workplace Cells.
- The installation and developer's guide for the i2b2 Workbench.
- The start up guide for the web client.
- The data and hive installation guides.
- The VM Install document.
- The CRC Plug-in Analysis document.
- A process document for building queries from ontology terms.

Fixes:

The following fixes were added to this release.

Query response messages remaining in queue

When a query "times out" a response message is filed in the JMS queue. These response messages are never cleared out of the queue. This can result in apparent prolonged response times.

STATUS/RESOLUTION: *This issue has been resolved. Query response messages will now be removed from the queue.*

The following outlines how to fix this issue in an existing system.

Edit File: JBOSS_HOME/server/default/deploy/crc-queryresponse-q-service.xml

Update lines:

```
<attribute name="RedeliveryLimit">1</attribute>
<attribute name="RedeliveryDelay">1</attribute>
```

Previous Query: Queries created in webclient will not run in workbench

Previous queries that were originally created in the webclient will not run in the workbench. When dragging the previous query to the Query Tool in the workbench an error would be returned saying you could not use the query.

STATUS/RESOLUTION: *This issue has been resolved. Queries created in either the webclient or workbench will run in both environments.*

The <observer_cd> in the PDO xml is missing the name attribute

In the PDO xml the xml tag for observer_cd is missing the name attribute. This is inconsistent with how the workbench works.

STATUS/RESOLUTION: *This issue has been resolved. The name attribute is now included.*

Licensing

The i2b2 source code is licensed under the i2b2 Software License 2.1. This includes but is not limited to all code in the edu.harvard.i2b2.* package namespace.

Trademarks

Microsoft, Windows, Windows XP, Windows Vista, Windows 7 and Windows Server are either registered trademarks or trademarks of Microsoft Corporation in the United States and/or other countries.

Mac OS is a registered trademark of Apple, Inc.

Eclipse Galileo is a registered trademark of the Eclipse Foundation, Inc.

All other trademarks are the property of their respective owners.