



343. Import DB Table

To import data from another sql database, you need to add the server as a Datasource Server first. Select Show Datasource Servers in the [Import Browser](#) and right click the free area and select Add Datasource Server.

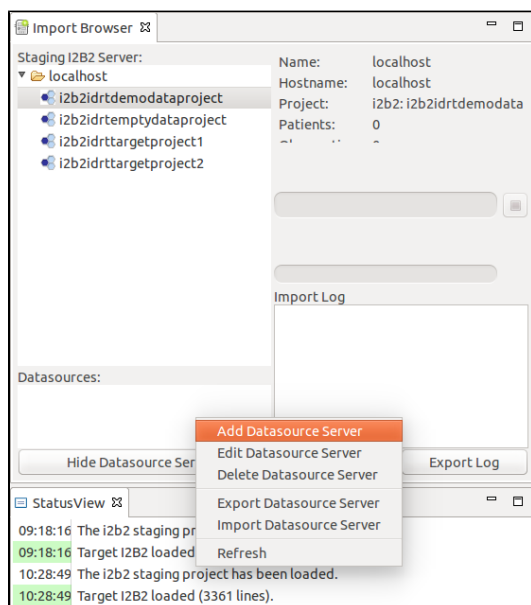

You can add a Datasource Server the same way, you added the [Staging Server](#).

To start the DB Table Import, right click the project in the [Import Browser](#) you want to upload data to and select Import Data -> Import DB Table.

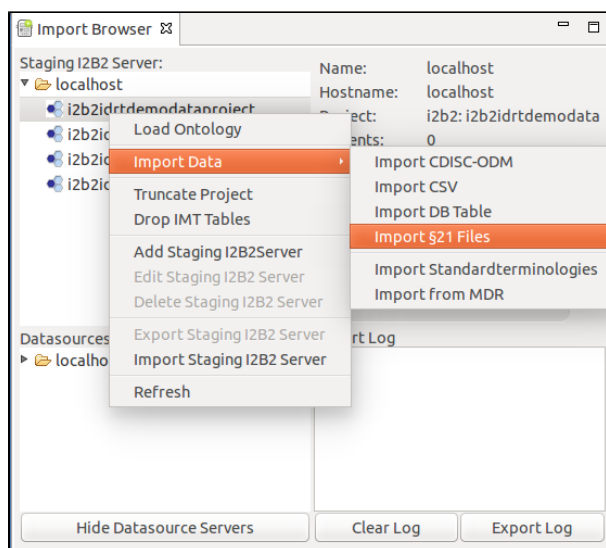

A window opens, where you can select the DB Table you want to import.

Source DB Settings

Error in table: db_import. PatientID missing?

☐ archive_observation_1
 ☐ birn
 ☐ code_lookup
 ☐ concept_dimension
 ☐ custom_meta
 ☐ datamart_report
 ☒ db_import
 ☐ encounter_mapping

Item	Name	Tooltip	Datatype	Metadata
PatientID			String	
Medication			String	
ICD10			String	

Guess Schema

Clear Single Table

Clear All Tables

< Back

Next >

Cancel

Finish

Item: The column name in the DB table.

Name: The name, the item should have in i2b2. You can e.g. rename cryptic column names.

Tooltip: The tooltip the item should have in i2b2. You can set a tooltip for better understanding, or the unit of an item for example.

Datatype: The datatype the item has. You can select either Integer, Float, Date, or String.

Metadata: This tells the Import Browser, where to look for specific items, such as the id of the patient, the import date, start date and so on. You can also ignore items. These will not end up in the i2b2 database.



The metadata item **PatientID/ObjectID** has to be present in each configuration file. Without it, the Import Browser won't let you continue and displays an error.



The **ObjectID** is the same as the **PatientID**, but the file is imported using modifiers instead. This way more than one row per patient is allowed. The Import Browser will tell you, if you used the Guess Schema functionality whether you should use the **ObjectID** instead.

The next table describes the buttons in the DB Import Settings page.

Name of the Button	Description
Clear Table	This will delete your configuration for the selected table.
Guess Schema	You can let the Import Browser try to guess the schema of the table. This will override your config!

If you select **Next** you get to the last page of the DB Importer. Here you can see all tables you selected and set some further options.

DB Import Settings

DB Import Settings

Truncate i2b2 Project? ☐

Truncate Previous Queries? ☐

Database Indexing ☒ Ignore ☐ Stop/Start ☐ Drop/Create

Set patient_count after import ☐

Save Settings? ☐

Tables to Import:

DS localhost: i2b2idrttargetproject2 - db_import

< Back Next > Cancel Finish

Name	Description
Truncate i2b2 Project?	If you check this, the project will be truncated before the data is uploaded.
Truncate Previous Queries?	This will truncate the previous queries tables.
Database Indexing	Ignore: Don't touch the indexes; Stop/Start: Stop the indexes before the upload and start them afterwards; Drop/Create: Drop the indexes before the upload and create them from scratch afterwards.
Set patient_count after import	This will fill the c_totalnum column.
Save Settings?	This will save the settings for the next import. Note: The truncate checks are never saved.

If you hit **Finish** the upload starts. You can observe the progress in the progress bar of the [Import Browser](#).